


## ADVANCED PHOTON SOURCE



# CAT Guide and Beamline Directory

A Key to APS Collaborative Access Teams

November 2000  
Argonne National Laboratory

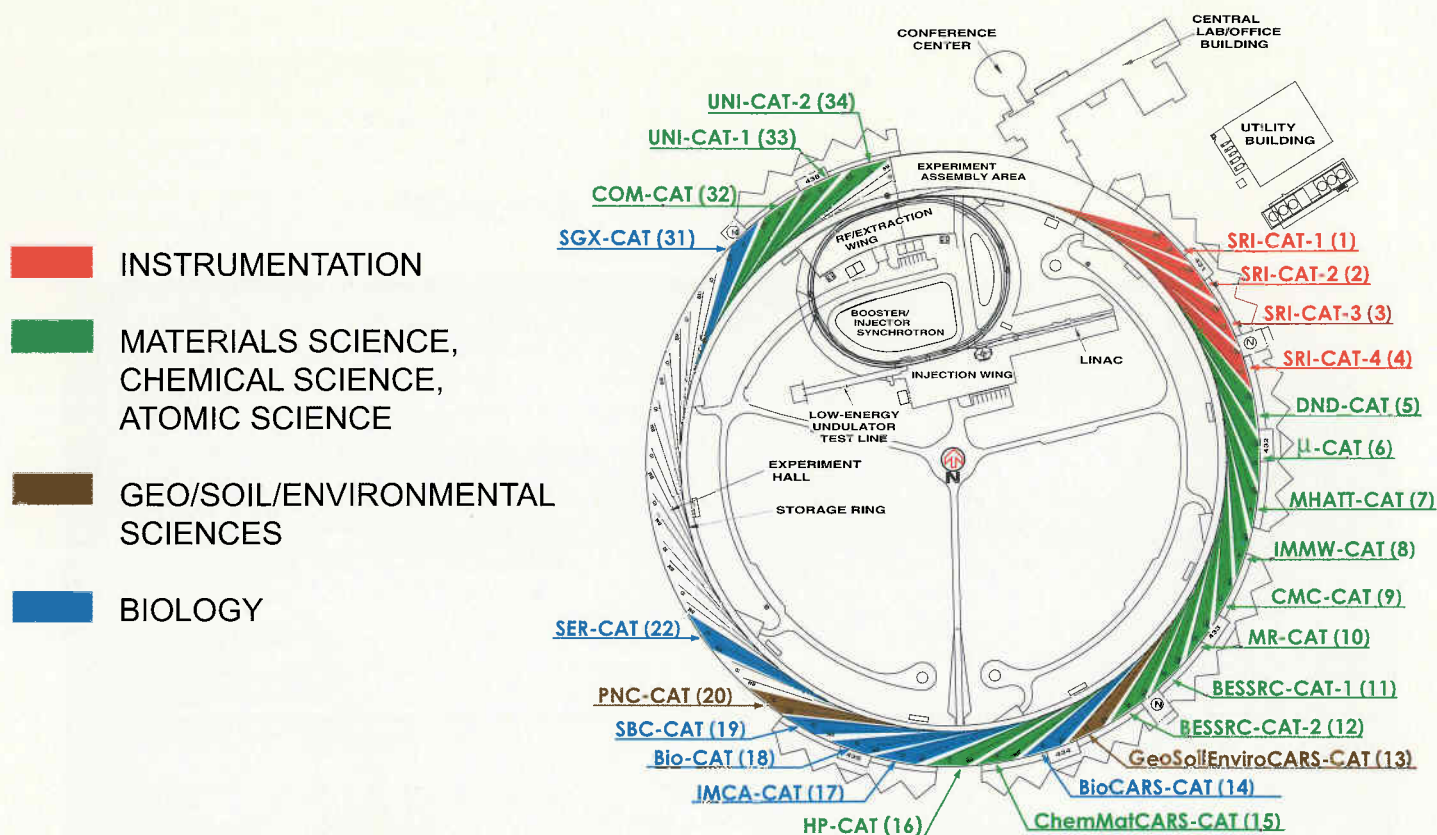

CAT ACRONYM	CAT NAME	SECTOR	BEAMLINE(S)
SRI-CAT	Synchrotron Radiation Instrumentation	1	1-BM, 1-ID
		2	2-BM, 2-ID
		3	3-ID
		4	4-ID
DND-CAT	DuPont–Northwestern–Dow	5	5-BM, 5-ID
μ-CAT	Midwest Universities	6	6-ID
MHATT-CAT	Center for Real-Time X-ray Studies	7	7-BM, 7-ID
IMMW-CAT	IBM–MIT–McGill–Whitehead	8	8-BM, 8-ID
CMC-CAT	Complex Materials Consortium	9	9-BM, 9-ID
MR-CAT	Materials Research	10	10-ID
BESSRC-CAT	Basic Energy Sciences Synchrotron Radiation Center	11	11-ID
		12	12-BM, 12-ID
CARS-CAT	Consortium for Advanced Radiation Sources		
	GeoSoilEnviroCARS	13	13-BM, 13-ID
	BioCARS	14	14-BM, 14-ID
	ChemMatCARS	15	15-ID
HP-CAT	High-Pressure	16	under development
IMCA-CAT	Industrial Macromolecular Crystallography Association	17	17-BM, 17-ID
Bio-CAT	Biophysics	18	18-BM, 18-ID
SBC-CAT	Structural Biology Center	19	19-BM, 19-ID
SER-CAT	Southeast Regional	22	under development
PNC-CAT	Pacific Northwest Consortium		
SGX-CAT	Structural GenomiX	31	under development
COM-CAT	Commercial	32	32-ID
UNI-CAT	A University–National Laboratory–Industry	33	33-BM, 33-ID
		34	34-ID

# CAT Guide and Beamline Directory

## TABLE OF CONTENTS

Facility Overview .....	1
User Program .....	1
Summary Beamline/Scientific Programs Table .....	2
SRI-CAT (Sectors 1, 2, 3, and 4) .....	7
DND-CAT (Sector 5) .....	17
$\mu$ -CAT (Sector 6) .....	21
MHATT-CAT (Sector 7) .....	23
IMMW-CAT (Sector 8) .....	25
CMC-CAT (Sector 9) .....	27
MR-CAT (Sector 10) .....	31
BESSRC-CAT (Sectors 11 and 12) .....	33
CARS-CAT (Sectors 13, 14, and 15) .....	37
HP-CAT (Sector 16) .....	45
IMCA-CAT (Sector 17) .....	49
Bio-CAT (Sector 18) .....	53
SBC-CAT (Sector 19) .....	55
PNC-CAT (Sector 20) .....	59
SER-CAT (Sector 22) .....	63
SGX-CAT (Sector 31) .....	67
COM-CAT (Sector 32) .....	69
UNI-CAT (Sectors 33 and 34) .....	71



Argonne National Laboratory is operated by The University of Chicago for the  
U.S. Department of Energy under Contract W-31-109-Eng-38.





# The Advanced Photon Source

---

## Facility Overview

The Advanced Photon Source (APS), a national user facility for synchrotron radiation research, is located at Argonne National Laboratory, approximately 25 miles southwest of Chicago, Illinois. The APS is considered a third-generation synchrotron radiation facility (specifically designed to accommodate insertion devices to serve as radiation sources) and is one of three such facilities in the world. Currently, it is the most brilliant source in the United States for research in such diverse fields as biology, medicine, materials science, chemistry, geology, agriculture and soil science, physics, and manufacturing technology.

Researchers use the APS either as members of Collaborative Access Teams (CATs) or as Independent Investigators (IIs). CATs are responsible for designing, building, and operating beamlines in one or more sectors, each sector consisting of an insertion-device (ID) beamline and a bending-magnet (BM) beamline. Each beamline is designed to accommodate a specific type of research program(s) and is optimized accordingly.

## User Program

CAT members are entitled to use 75% of the available beam time to pursue CAT research goals. The remaining 25% of the available beam time must be made available to IIs, who are selected by the CAT through a proposal process determined by each CAT and approved by the APS. Prospective IIs submit proposals to the APS User Office via the World Wide Web (WWW) at the following address: [http://www.aps.anl.gov/xfd/communicator/useroffice/II\\_proposal.html](http://www.aps.anl.gov/xfd/communicator/useroffice/II_proposal.html). Each proposal is electronically logged in and reviewed according to the review process designated by the selected CAT. This document, the *CAT Guide and Beamline Directory*, was written to help prospective IIs determine which beamlines are suitable for their specific experiments. Further information can be obtained from the APS User Office or the individual CATs (contact information is provided on the inside back cover of this document).

# Summary Beamline / Scientific Programs Table

BEAMLINE	CAT	SCIENTIFIC FOCUS	PROGRAMS	SOURCE	STATUS *
1-BM	SRI	Synchrotron instrumentation and techniques, time-resolved studies	-Time-resolved studies, diffraction, and spectroscopy	BM	O II
1-ID	SRI	Synchrotron instrumentation and techniques, high-energy x-rays, optics development	-High-energy x-ray scattering techniques -High-heat-load x-ray optics development	Und A (3.3 cm)	O II
2-BM	SRI	Synchrotron instrumentation and techniques, x-ray lithography, microtomography	-Deep x-ray lithography -X-ray microtomography -Development of x-ray optics -Beam diagnostics instrumentation -X-ray microprobe studies	BM	O II
2-ID-B	SRI	Synchrotron instrumentation and techniques, imaging, coherence	-High-resolution imaging -Coherent scattering -Coherence	5.5-cm Und	O II
2-ID-D/E	SRI	Synchrotron instrumentation and techniques, imaging, coherence	-High-resolution imaging -Coherence-based techniques	Und A (3.3 cm)	O II
3-ID	SRI	Synchrotron instrumentation and techniques, resonant and non-resonant inelastic scattering	-Nuclear resonant scattering -Inelastic x-ray scattering	2.7-cm Und	O II
4-ID	SRI	Synchrotron instrumentation and techniques, production and use of polarized x-rays	-Development of polarization optics and techniques for 0.5 - 100 keV -Development of high-heat-load front-end components	Und A (3.3 cm) & Circularly Polarized Und	D
5-BM	DND	Crystallography, EXAFS, polymer scattering, tomography	-Materials science and engineering -Polymer science -Crystallography -Surface science -Catalysis -Imaging	BM	O C II
5-ID	DND	Crystallography, EXAFS, polymer scattering, tomography	-Materials science and engineering -Polymer science -Crystallography -Surface science -Catalysis	Und A (3.3 cm)	O C II

6-ID	$\mu$	Materials science	-Magnetic scattering and spectroscopy -Liquid-surface studies -Surface scattering	Und A (3.3 cm)	C
7-BM	MHATT	Real-time x-ray studies	-Real-time x-ray studies -X-ray microcharacterization -X-ray photon correlation spectroscopy	BM	D
7-ID	MHATT	Real-time x-ray studies	-Real-time x-ray studies -X-ray microcharacterization -X-ray photon correlation spectroscopy	Und A (3.3 cm)	O C
8-ID	IMMW	Condensed matter physics, materials science	-Intensity fluctuation spectroscopy studies using coherent x-rays -Time-resolved x-ray scattering -X-ray scattering studies at very-low temperatures -X-ray magnetic scattering	Und A (3.3 cm)	O C
9-BM	CMC	Materials science	-EXAFS of complex materials -X-ray surface scattering -General diffraction	BM	D
9-ID	CMC	Materials science	-SAXS from complex materials -Liquid/solid surface scattering -EXAFS -General diffraction -Magnetic and inelastic	Und A (3.3 cm)	C
10-ID	MR	Materials science, environmental science	-Spectroscopy of dilute environmental systems -Combined techniques for <i>in situ</i> characterization of materials -Zone plate microfocusing spectroscopy and imaging	Und A (3.3 cm)	O C
11-ID	BESSRC	Spectroscopy, scattering, high-energy diffraction, magnetic Compton scattering	-Materials science: x-ray scattering, high-energy diffraction, dichroism, magnetism -Geoscience -Atomic physics -Chemistry: spectroscopy (EXAFS, XANES), time-dependent spectroscopy of excited states	Elliptical Multipole Wiggler	O C II

\* O = operational (beamlines identified as operational may have facilities still being commissioned.)

C = commissioning

D = design / construction

II = accepting Independent Investigator proposals

BM = bending magnet

Und = undulator

BEAMLINE	CAT	SCIENTIFIC FOCUS	PROGRAMS	SOURCE	STATUS *
12-BM	BESSRC	Spectroscopy, x-ray scattering, diffraction	<ul style="list-style-type: none"> <li>-Materials science: x-ray scattering, diffraction and spectroscopy</li> <li>-Geoscience: mineral-liquid interface characterization by x-ray scattering and spectroscopy</li> <li>-Chemistry: EXAFS and XANES spectroscopy of actinides and lanthanides</li> <li>-Atomic physics</li> </ul>	BM	O C II
12-ID	BESSRC	Materials science, atomic physics, geoscience	<ul style="list-style-type: none"> <li>-Materials science: elastic and inelastic scattering, MOCVD, x-ray standing waves, electrochemistry</li> <li>-Geoscience: mineral-liquid interface characterization</li> <li>-Chemistry: small-angle scattering</li> <li>-Atomic physics, spectroscopy</li> </ul>	Und A (3.3 cm)	O C II
13-BM	CARS GSECARS	Geosciences, environmental science, soil science	<ul style="list-style-type: none"> <li>-High-pressure diffraction in diamond-anvil cell</li> <li>-High-pressure diffraction in multi-anvil press</li> <li>-Microspectroscopy and XRF microprobe</li> <li>-Microtomography</li> </ul>	BM	O
13-ID	CARS GSECARS	Geosciences, environmental science, soil science	<ul style="list-style-type: none"> <li>-High-pressure diffraction in diamond-anvil cell</li> <li>-High-pressure diffraction in multi-anvil press</li> <li>-Microspectroscopy and XRF microprobe</li> <li>-Microcrystal diffraction</li> <li>-Surface diffraction</li> <li>-X-ray absorption fine structure spectroscopy</li> </ul>	Und A (3.3 cm)	O
14-BM	CARS BioCARS	Structural biology	<ul style="list-style-type: none"> <li>-Large unit cell (virus) crystallography</li> <li>-Small unit cell (protein) crystallography</li> <li>-MAD phasing</li> <li>-Time-resolved crystallography</li> <li>-Laue diffraction</li> <li>-Study of microcrystals</li> </ul>	BM	O



14-ID	CARS BioCARS	Structural biology	-See 14-BM	Wiggler A (8.5 cm)	O
15-ID	CARS ChemMat- CARS	Dynamic and structural condensed matter chemistry	-Time-dependent chemical crystallography scattering -Static and time-dependent surface scattering studies of surface and interfacial properties -Small-angle and wide-angle x-ray scattering -Coherent x-ray scattering studies of dynamic properties	Und A (3.3 cm)	C
16-BM	HP	High-pressure research	-White beam and monochromatic diffraction -EXAFS -Techniques development	BM	D
16-ID	HP	High-pressure research	-Angle-dispersive microdiffraction -Nuclear resonance spectroscopy -Electronic spectroscopies -Single-crystal diffraction -Inelastic scattering -Techniques development	Und A (3.3 cm)	D
17-BM	IMCA	Pharmaceutical macromolecular crystallography	-Structures of protein-ligand complexes - <i>De novo</i> protein structures -Drug design -Protein engineering -Crystallographic methods development	BM	O
17-ID	IMCA	Pharmaceutical macromolecular crystallography	-See above	Und A (3.3 cm)	O II
18-ID	Bio	Static and dynamic studies of partially ordered biological systems	-Fiber diffraction -X-ray scattering -X-ray absorption/emission spectroscopy -X-ray imaging and tomography	Und A (3.3 cm)	O
19-BM	SBC	Macromolecular crystallography	-Structural biology	BM	O
19-ID	SBC	Macromolecular crystallography	-Structural biology	Und A (3.3 cm)	O

\* O = operational (beamlines identified as operational may have facilities still being commissioned.)

C = commissioning

D = design / construction

II = accepting Independent Investigator proposals

BM = bending magnet

Und = undulator

BEAMLINE	CAT	SCIENTIFIC FOCUS	PROGRAMS	SOURCE	STATUS *
20-BM	PNC	Materials science, environmental science	-Materials science -Environmental science -Surface science	BM	C
20-ID	PNC	Materials science, environmental science	-Materials science -Environmental science -UHV/surface science, MBE growth -X-ray microbeams	Und A (3.3 cm)	O
22-BM	SER	Macromolecular crystallography	-Structural biology	BM	D
22-ID	SER	Macromolecular crystallography	-Structural biology	Und A (3.3 cm)	D
31-BM	SGX	Structural genomics	-Under development	BM	D
31-ID	SGX	Structural genomics	-Under development	Und A (3.3 cm)	D
32-ID	COM	Materials characterization	-Analytical materials characterization -"Fee-for-service" analyses	Und A (3.3 cm)	D
33-BM	UNI	Materials science and condensed matter physics	-Materials science -Ceramic science -Phase transitions -Surface science -Thin-films structure and growth -Materials physics	BM	D
33-ID	UNI	Materials science and condensed matter physics	-Materials science -Ceramic science -Phase transitions -Surface science -Thin-films structure and growth -Materials physics	Und A (3.3 cm)	O II
34-ID	UNI	Microfocusing and coherent x-ray scattering	-Coherent x-ray diffraction -Microbeam scattering	Und A (3.3 cm)	D

\* O = operational (beamlines identified as operational may have facilities still being commissioned.)

C = commissioning

D = design / construction

II = accepting Independent Investigator proposals

BM = bending magnet

Und = undulator



## Synchrotron Radiation Instrumentation CAT

### Overview

SRI-CAT is organized into developers and scientific members, with the developers coming from the Experimental Facilities Division and the User Program Division of the APS, as well as from Purdue University, University of Houston, National Institute of Science and Technology, Lawrence Berkeley National Laboratory, Stanford Synchrotron Radiation Laboratory, Cornell High Energy Synchrotron Source, and many Australian institutions managed by the Australian Nuclear Science and Technology Organization. Scientific members whose interests match those of the developers are drawn from the synchrotron user community.

### Research Focus

The principal mission of SRI-CAT is to develop new and unique forefront instruments and techniques to advance the use of synchrotron radiation. To that end, research programs involve the design, construction, and operation of standard and specialized insertion devices; diagnostics of radiation from these insertion devices; development of high-heat-load optics; the design, construction, and testing of various standard and modular beamline components; and the development of strategic instruments and techniques.

Among those instruments/capabilities are microbeam techniques, nuclear resonant spectroscopy, coherence-based techniques, the development of 0.5 to 4 keV scientific capabilities, the development of scientific capabilities at high energies ( $> 50$  keV), development of time-resolved instrumentation and techniques, the development of ultrahigh-resolution inelastic scattering techniques, and the development of instrumentation and techniques for application of variably polarized x-rays from 0.5 keV to 100 keV.

---

#### CAT contacts:

Dennis Mills, Executive Director  
Efim Gluskin, Program Director  
Trudy Bolin, CAT Manager  
Linda Shoudis, Admin. Asst.

tel 630.252.5680  
tel 630.252.4788  
tel 630.252.0801  
tel 630.252.0160

dmm@aps.anl.gov  
gluskin@aps.anl.gov  
bolitru@aps.anl.gov  
ljs@aps.anl.gov

# Beamline 1-BM / SRI-CAT

**Scientific focus:** Synchrotron instrumentation and techniques and time-resolved studies

**Scientific programs:** Time-resolved studies, dispersive diffraction, and spectroscopy

## Bending magnet Source Characteristics (nominal)

source	APS bending magnet
critical energy	19.51 keV
energy range	1.0 - 100.0 keV
on-axis peak brilliance at 16.3 keV	$2.9 \times 10^{15}$ ph/sec/mrad <sup>2</sup> /mm <sup>2</sup> /0.1%bw
on-axis peak angular flux at 16.3 keV	$9.6 \times 10^{13}$ ph/sec/mrad <sup>2</sup> /0.1%bw
on-axis peak horizontal angular flux at 5.6 keV	$1.6 \times 10^{13}$ ph/sec/mradh/0.1%bw
source size at critical energy $\sum x$ $\sum y$	145 $\mu$ m 36 $\mu$ m
source divergence at critical energy $\sum x'$ $\sum y'$	6 mrad 47 $\mu$ rad
1-BM acceptance	3.7 mrad

## Optics & Optical Performance

- standard mode
  - vertical collimating mirror
  - 25.5 m from source
  - bent cylinder
- PSL double-crystal monochromator
  - 27.2 m from source
  - sagittal focus
  - Si(111) crystal
  - $1.5 \times 10^{-4}$  energy resolution ( $\Delta E/E$ )
  - 35 mm beam offset
  - water cooling
  - focusing  $\approx 0.25 \times 0.6$  mm
  - $9 \times 10^{11}$  (photons/sec)/100mA at 10keV
- vertical focusing mirror
  - 45.5 m from source
  - cylindrical bend
- dispersive crystal mode
  - bent-crystal monochromator
  - 31.5 m from source
  - horizontal focus
  - with or without first mirror

## Experiment Stations

- 1-BM-A
  - first optics enclosure
  - collimating mirror
  - PSL monochromator
- 1-BM-B
  - white, pink, and monochromatic beam station
  - dispersive modes
- 1-BM-C
  - monochromatic beam station

## Detectors

- NaI scintillation counters
- ionization chambers
- Si(Li) and Ge energy dispersive detectors
- CCDs
- image plates

## Beamline Controls and Data Acquisition

- Sun UNIX running EPICS with VME
- SPEC

## Beamline Support Equipment/Facilities

- 6-circle Huber diffractometer
- image plate reader

**Beamline contacts:** Dean Haeffner  
Peter Lee  
Jin Wang

tel 630.252.0126  
tel 630.252.0162  
tel 630.252.9125

haeffner@aps.anl.gov  
pllee@aps.anl.gov  
wang@aps.anl.gov

# Beamline 1-ID / SRI-CAT

**Scientific focus:** Synchrotron instrumentation and techniques, high-energy x-rays, and optics development

**Scientific programs:** High-energy x-ray scattering techniques and x-ray optics development

## Insertion-device Source Characteristics (nominal)

source	Undulator A
period	3.30 cm
length	2.47 m
effective $K_{\max}$ (at minimum gap = 10.5 mm)	2.78
energy range 1st harmonic	2.9 - 13.0 keV
energy range 1st - 5th harmonics	2.9 - 45.0 keV
on-axis peak brilliance at 6.5 keV	$9.6 \times 10^{18}$ ph/sec/mrad <sup>2</sup> /mm <sup>2</sup> /0.1% bw
source size at 8.0 keV $\sum_x$ $\sum_y$	359 $\mu\text{m}$ 21 $\mu\text{m}$
source divergence at 8.0 keV $\sum_{x'}$ $\sum_{y'}$	24 $\mu\text{rad}$ 6.9 $\mu\text{rad}$

## Optics & Optical Performance

- Kohzu double-crystal constant off-set monochromator
  - 30 m from source
  - 35 mm offset
  - liquid-nitrogen cooling
  - 6–19 keV energy range Si(111)
  - $1.8 \times 10^{13}$  photons/sec at 10 keV on sample (unfocused)
  - $10^{-4}$  energy resolution ( $\Delta E/E$ ) at 10 keV
  - 10–42 keV energy range Si(311)
  - $7 \times 10^{11}$  photons/sec at 19.2 keV,  $\Delta E = 2.5$  eV
  - $3 \times 10^{11}$  photons/sec at 33.2 keV,  $\Delta E = 4.5$  eV
- high-energy monochromator
  - 35 m from source
  - Si(111) limited tunability
  - 50–150 keV energy range
  - $1 \times 10^{12}$  photons/sec at 100 mA (1 x 1 mm<sup>2</sup> at 60 m)
  - $\Delta E/E = 1 \times 10^{-3}$

## Experiment Stations

- 1-ID-A
  - white beam first optics enclosure
  - Kohzu monochromator
- 1-ID-B
  - white beam station
  - high-energy monochromator
- 1-ID-C
  - monochromatic beam station

## Detectors

- NaI scintillation
- ionization chambers
- Princeton Instruments CCD
- Ge & Si(Li) solid state detectors
- image plate

## Beamline Controls and Data Acquisition

- Sun UNIX running EPICS with VME
- SPEC

## Beamline Support Equipment/Facilities

- 6-circle diffractometer

**Beamline contacts:** Dean Haeffner  
Wah Keat Lee  
Peter Lee  
Sarvjit Shastri

tel 630.252.0126  
tel 630.252.7759  
tel 630.252.0162  
tel 630.252.0129

haeffner@aps.anl.gov  
wklee@aps.anl.gov  
pllee@aps.anl.gov  
shastri@aps.anl.gov

# Beamline 2-BM / SRI-CAT

**Scientific focus:** Synchrotron instrumentation and techniques, x-ray lithography, and microtomography

**Scientific programs:** Deep x-ray lithography, x-ray microtomography, development of x-ray optics, beam diagnostics instrumentation, and x-ray microprobe studies

## Bending Magnet Source Characteristics (nominal)

source	APS bending magnet
critical energy	19.51 keV
energy range	1.0 - 100.0 keV
on-axis peak brilliance at 16.3 keV	$2.9 \times 10^{15}$ ph/sec/mrad <sup>2</sup> /mm <sup>2</sup> /0.1%bw
on-axis peak angular flux at 16.3 keV	$9.6 \times 10^{13}$ ph/sec/mrad <sup>2</sup> /0.1%bw
on-axis peak horizontal angular flux at 5.6 keV	$1.6 \times 10^{13}$ ph/sec/mradh/0.1%bw
source size at critical energy $\sum_x$ $\sum_y$	145 $\mu$ m 36 $\mu$ m
source divergence at critical energy $\sum_{x'}$ $\sum_{y'}$	6 mrad 47 $\mu$ rad
2-BM acceptance	2.0 mrad

## Optics & Optical Performance

- vertical deflecting mirror  
24.9 m from source  
0.15° incident angle  
Cr, Pt coatings  
plane figure
- double multilayer monochromator  
27.4 m from source  
unfocussed
- Kohzu double-crystal monochromator  
28.8 m from source  
unfocussed
- vertical deflecting mirror  
49.4 m from source  
0°–2° incident angle, variable  
Pt coating  
plane figure

## Experiment Stations

- 2-BM-A  
white beam first optics enclosure
- 2-BM-B  
pink and monochromatic beam station

## Detectors

- energy dispersive detector
- scintillation detector
- ionization chambers
- peltier-cooled CCD camera

## Beamline Controls and Data Acquisition

- Sun UNIX running EPICS with VME; NT and MacOS workstations
- SPEC

## Beamline Support Equipment/Facilities

- precision scanning stage for x-ray lithography
- 6-circle diffractometer
- first x-ray tomography system

**Beamline contacts:** Derrick Mancini  
Francesco DeCarlo

tel 630.252.0147  
tel 630.252.0148

mancini@aps.anl.gov  
decarlo@aps.anl.gov



# Beamline 2-ID-B / SRI-CAT

**Scientific focus:** Synchrotron instrumentation and techniques, imaging, and coherence

**Scientific programs:** High-resolution imaging, coherent scattering, and interferometry

## Insertion Device Source Characteristics (nominal)

source	5.5 cm undulator
period	5.50 cm
length	2.47 m
effective $K_{\max}$ (at minimum gap = 10.5 mm)	6.57
energy range 1st harmonic	0.4 - 7.0 keV
energy range 1st - 5th harmonics	0.4 - 25.0 keV
on-axis peak brilliance at 4.0 keV	$4.2 \times 10^{18}$ ph/sec/mrad <sup>2</sup> /mm <sup>2</sup> /0.1%bw
source size at 4.0 keV $\sum_x$ $\sum_y$	359 $\mu\text{m}$ 21 $\mu\text{m}$
source divergence at 4.0 keV $\sum_{x'}$ $\sum_{y'}$	24 $\mu\text{rad}$ 9.0 $\mu\text{rad}$

## Optics & Optical Performance

- high-heat-load mirror M1
  - 29.5 m from source
  - 0.15° incident angle
  - plane figure
  - Pt, Rh, and Si coatings
- horizontal focusing mirror M2B
  - 31.1 m from source
  - 1.25° incident angle
  - spherical figure
  - Pt, Rh, multilayer coatings
- vertical focusing mirror M3B
  - 31.6 m from source
  - 1.25° incident angle
  - cylindrical figure
  - Pt, Rh, multilayer coatings
- multilayer spherical grating monochromator
  - 46.8 m from source
  - 0.5–4.0 keV energy range
  - 600, 1200, and 1800 lines/mm gratings
  - 100–4000 resolution  $E/\Delta E$
  - 1.5 mm hor. x 0.5 mm vert. beam size FWHM

- microprobe zone plate
  - 60.3 m from source
  - 60 nm x 60 nm focus size
  - $10^7$ – $10^8$  ph/sec/0.1% bw microprobe flux

## Experiment Station

- 2-ID-B

## Detectors

- absolute-calibrated photodiodes
- avalanche photodiodes
- dispersive LEGe fluorescence detector
- thinned/backside-illuminated CCD camera

## Beamline Controls and Data Acquisition

- Sun UNIX running EPICS with VME
- PC NT running EPICS CCD camera controller

## Beamline Support Equipment/Facilities

- scanning/holographic microscope
- 2-circle goniometer
- fast scan stage (0.8 nm resolution)
- full-field imaging microscope

**Beamline contacts:** Ian McNulty  
Sean Frigo

tel 630.252.2882  
tel 630.252.9374

mcnulty@aps.anl.gov  
frigo@aps.anl.gov

# Beamline 2-ID-D & 2-ID-E / SRI-CAT

**Scientific focus:** Synchrotron instrumentation and techniques, imaging, and coherence

**Scientific programs:** High-resolution imaging and coherence-based techniques

## Insertion-device Source Characteristics (nominal)

source	Undulator A
period	3.30 cm
length	2.47 m
effective $K_{\max}$ (at minimum gap = 10.5 mm)	2.78
energy range 1st harmonic	2.9 - 13.0 keV
energy range 1st - 5th harmonics	2.9 - 45.0 keV
on-axis peak brilliance at 6.5 keV	$9.6 \times 10^{18}$ ph/sec/mrad <sup>2</sup> /mm <sup>2</sup> /0.1% bw
source size at 8.0 keV $\sum x$ $\sum y$	359 $\mu\text{m}$ 21 $\mu\text{m}$
source divergence at 8.0 keV $\sum x'$ $\sum y'$	24 $\mu\text{rad}$ 6.9 $\mu\text{rad}$

## Optics & Optical Performance

- high-heat-load mirror M1  
29.5 m from source  
0.15° incident angle  
plane figure  
Pt, Rh, and Si coatings
- Kohzu monochromator  
62.0 m from source  
2–32 keV energy range Si(111)  
10<sup>-4</sup> monochromaticity  $\Delta E/E$   
4.2 mm hor. x 1.6 mm vert. beam size FWHM  
10<sup>11</sup>–10<sup>13</sup> ph/sec flux at sample
- zone plate microprobe  
71.0 m from source  
0.1  $\mu\text{m}$  hor. x 0.1  $\mu\text{m}$  vert. focus size  
10<sup>11</sup> ph/sec/ $\mu\text{m}^2$ /0.1% bw flux

## Experiment Stations

- 2-ID-D
- 2-ID-E

## Detectors

- ionization chambers
- NaI scintillation detector
- avalanche photodiodes
- CCD cameras

## Beamline Controls and Data Acquisition

- Sun UNIX running EPICS with VME
- PC NT running EPICS CCD camera controller

## Beamline Support Equipment/Facilities

- optical table with six degrees of freedom
- scanning microscope
- high-resolution x-ray spectrometer
- Kappa diffractometer

## Beamline contacts:

Barry Lai  
Zhonghou Cai  
Jorg Maser

tel 630.252.6405  
tel 630.252.0144  
tel 630.252.1081

blai@aps.anl.gov  
cai@aps.anl.gov  
maser@aps.anl.gov

# Beamline 3-ID / SRI-CAT

**Scientific focus:** Synchrotron instrumentation and techniques, resonant and non-resonant inelastic scattering

**Scientific programs:** Nuclear resonant scattering (inelastic nuclear resonant scattering, lattice dynamics of thin films, amorphous materials, and vibrational properties of biomolecules), and inelastic x-ray scattering (momentum-resolved lattice dynamic)

## Insertion Device Source Characteristics (nominal)

source	2.7 cm undulator
period	2.70 cm
length	4.5 m
effective $K_{\text{max}}$ (at minimum gap = 8.5 mm)	2.18
energy range 1st harmonic	5.1 - 16.0 keV
energy range 1st - 5th harmonics	5.1 - 60.0 keV
on-axis peak brilliance at 7.9 keV	$2.3 \times 10^{19}$ ph/sec/mrad <sup>2</sup> /mm <sup>2</sup> /0.1%bw
source size at 8.0 keV $\sum_x$ $\sum_y$	359 $\mu\text{m}$ 21 $\mu\text{m}$
source divergence at 8.0 keV $\sum_{x'}$ $\sum_{y'}$	24 $\mu\text{rad}$ 6.9 $\mu\text{rad}$

## Optics & Optical Performance

- Kohzu HLM-2 double-crystal monochromator  
diamond (111) nondispersive crystals  
water cooling
- high-resolution monochromator
 

energy	resolution
13.8 keV	6 meV
14.4 keV	5.5, 2, 0.85, and 0.66 meV
21.5 keV	1 meV
23.8 keV	3.5, 1, and 0.5 meV
- Zeiss mirror  
Zerodur substrate  
Pd coating  
two groove, collimating and horizontal  
focusing  
bendable for vertical focusing  
0.8 m length
- K-B focusing mirror

## Experiment Stations

- 3-ID-A  
white beam first optics enclosure
- 3-ID-B  
monochromatic beam station  
x-ray optics development  
nuclear resonant scattering
- 3-ID-C  
monochromatic beam station  
inelastic x-ray scattering
- 3-ID-D  
monochromatic beam station  
nuclear resonant scattering

## Detectors

- avalanche photodiodes
- PIN-diodes, Si, CdZnTe (Amptek)
- ionization chambers
- NaI
- Si(Li)

## Beamline Controls and Data Acquisition

- VME-based motion control, EPICS software
- Graphical user interface: MEDM implemented  
on Sun workstations, Solaris operating system

## Beamline Support Equipment/Facilities

- Kohzu high-angular resolution stages (12.5 nrad)
- 6-circle diffractometer
- optical tables
- Euler cradle
- misc. linear stages
- Oxford 7-T superconducting magnet/cryostat,  
Oxford microscope cryostat
- optical tables

**Beamline contacts:** Ercan Alp  
Wolfgang Sturhahn

tel 630.252.4775  
tel 630.252.0163

eea@aps.anl.gov  
sturhahn@aps.anl.gov

# Beamline 4-ID / SRI-CAT

**Scientific focus:** Synchrotron instrumentation and techniques and use of polarized x-rays

**Scientific programs:** Development of polarization optics/techniques for 0.5–100 keV (magnetic circular dichroism, resonant magnetic scattering, and fluorescence), and development of high-heat-load front-end components.

## Circularly Polarized Undulator (nominal)\*

period	12.8 cm
length	2.3 m
effective $K_{\text{max}}$ (for both horizontal and vertical fields)	2.65
peak field $B_{\text{max}}$	0.28 Tesla
maximum currents at $B_{\text{max}}$	1.4 kA horizontal 0.32 kA vertical
energy range 1st harmonic (helical mode)	0.44 - 3.0 keV
energy range 1st harmonic (linear mode)	0.8 - 3.0 keV
energy range 1st - 5th harmonics (linear mode)	0.8 - 10.0 keV
on-axis peak circular brilliance at 1.5 keV	$1.0 \times 10^{18}$ ph/sec/mrad <sup>2</sup> /mm <sup>2</sup> /0.1%bw
on-axis peak linear brilliance at 1.9 keV	$7.0 \times 10^{17}$ ph/sec/mrad <sup>2</sup> /mm <sup>2</sup> /0.1%bw
switching frequency	0 - 10 Hz
switching rise time	20 ms
source size at 1.5 keV $\sum_x$ $\sum_y$	359 $\mu$ m 21 $\mu$ m
source divergence at 1.5 keV $\sum_{x'}$ $\sum_{y'}$	27 $\mu$ rad 14.7 $\mu$ rad

\* Sector 4 will be equipped with two canted insertion devices, a standard Undulator A and an elliptically polarized undulator.

## Optics & Optical Performance

Hard x-ray branch optics and optics performance

- double-crystal monochromator
  - 3.0–50.0 keV energy range Si(111)
  - 5 mm hor. x 2 mm vert. acceptance
  - 1°–60° Bragg angle rotation range
  - 10 arcsec Bragg accuracy
  - 0.5 arcsec resolution
  - 10–35 mm beam offset, variable
  - liquid-nitrogen cooling

## Insertion Device Source Characteristics (nominal)\*

source	Undulator A
period	3.30 cm
length	2.47 m
effective $K_{\text{max}}$ (at minimum gap = 10.5 mm)	2.78
energy range 1st harmonic	2.9 - 13.0 keV
energy range 1st - 5th harmonics	2.9 - 45.0 keV
on-axis peak brilliance at 6.5 keV	$9.6 \times 10^{18}$ ph/sec/mrad <sup>2</sup> /mm <sup>2</sup> /0.1% bw
source size at 8.0 keV $\sum_x$ $\sum_y$	359 $\mu$ m 21 $\mu$ m
source divergence at 8.0 keV $\sum_{x'}$ $\sum_{y'}$	24 $\mu$ rad 6.9 $\mu$ rad

### • mirror

two groove, torroidal figure  
80.5 mm sagittal radius, focusing  
345 mm sagittal radius, collimating  
18–7 km meridional radius, adjustable  
800 mm L x 30 mm W  
1.5–3.85 mrad angular range

Intermediate x-ray branch optics and optics performance

- horizontal focusing mirror M1C
  - 28.6 m from source
  - 1.1° incident angle
  - plane figure
  - Pt, Rh, Si, and multilayer coatings
- horizontal focusing mirror M2C
  - 31 m from source
  - 1.1° incident angle
  - spherical figure (R=1610 m)
  - Pt, Rh, Si, and multilayer coatings

*continued on pg. 15*

**Beamline contacts:** John Freeland  
George Srajer

tel 630.252.9614  
tel 630.252.3257

freeland@aps.anl.gov  
srajerg@aps.anl.gov

## Beamline 4-ID / SRI-CAT

**Scientific focus:** Synchrotron instrumentation and techniques and use of polarized x-rays

**Scientific programs:** Development of polarization optics/techniques for 0.5–100 keV (magnetic circular dichroism, resonant magnetic scattering, and fluorescence), and development of high-heat-load front-end components.

*continued from pg. 14*

- vertical focusing mirror M3C
  - 41.4 m from source
  - 1.00° incident angle
  - spherical figure
  - Rh coating
- multilayer spherical grating monochromator
  - 50 m from source
  - 0.5–3.0 keV energy range
  - $10^{-3}$ – $10^{-4}$  monochromaticity  $\Delta E/E$
  - 2.5 mm hor. x 0.25 mm vert. beam size
  - $10^{11}$ – $10^{13}$  ph/sec flux at sample

### Experiment Stations

- 4-ID-A
  - white beam first optic enclosure
- 4-ID-B and -D
  - white and monochromatic “hard” x-ray stations
- 4-ID-C
  - soft x-ray station

### Beamline Controls and Data Acquisition

- Sun UNIX running Epics with VME, SPEC


**Beamline contacts:** John Freeland  
George Srajer

tel 630.252.9614  
tel 630.252.3257

freeland@aps.anl.gov  
srajerg@aps.anl.gov





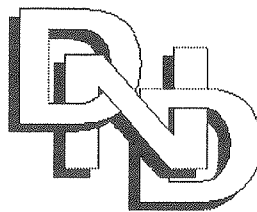

## **E.I. Du Pont de Nemours & Co-Northwestern University-The Dow Chemical Company CAT**

---

### **Overview**

The DND-CAT is a team of senior scientists and engineers formed jointly by the E.I. Du Pont de Nemours Co., Northwestern University, and the Dow Chemical Company to develop, build, and operate facilities dedicated to advancing x-ray research on new materials. The ability to characterize the increasing numbers of new and advanced materials with specifically engineered properties is closely coupled with the ability to analyze their molecular structure.

### **Research Focus**

The scientific thrust of DND-CAT is concentrated in three main areas: the study of the atomic structure of bulk materials, the study of two-dimensional atomic structures (e.g., surfaces, interfaces, and thin films), and polymer science and technology.

#### *Atomic Structure of Bulk Materials*

In addition to crystallographic determinations of unknown compounds, this aspect of DND-CAT's program concentrates on the determination of defects and local atomic arrangements in materials. For example, the arrangement of dopant atoms and other lattice defects in metals, catalysts, and semiconductors is under study. In addition, DND-CAT researchers are extensively studying catalysts under realistic conditions of high temperature and pressure. Catalysts of interest include supported metal and bimetallic catalysts (e.g., Pd, Pd-Re), zeolites, and catalysts for environmentally benign chlorofluorocarbons.

#### *Study of Two-Dimensional Structures*

The structure of surfaces at interfaces is a primary factor controlling the properties of materials in many technologically important applications. Some of these include catalysis, corrosion, tribology, crystal growth, and ion transport through biomembranes. DND-CAT researchers will use a custom-designed ultra-high-vacuum system to study buried interfaces, thin films, and surfaces to determine mechanisms of film growth, segregation at grain boundaries, surface melting, and other two-dimensional phenomena.

#### *Polymer Science and Engineering*

The third major area for DND-CAT is polymer physics, chemistry, and engineering. Areas of interest include solution characterization of polymers and their precursors, structure formation, polymer melt rheology, fiber formation, structure/property/processing relationships, and mechanical properties.

---

**CAT contacts:**

John Quintana, Director  
Jerrie Lea Hopf, Program Asst.  
WWW Site

tel 630.252.0221  
tel 630.252.0222

jpq@nwu.edu  
jlhopf@nwu.edu  
<http://tomato.dnd.aps.anl.gov/>

# Beamline 5-BM / DND-CAT

**Scientific focus:** Crystallography, EXAFS, polymer scattering, and tomography

**Scientific programs:** Material science and engineering, polymer science, crystallography, surface science, catalysis, and imaging

## Bending Magnet Source Characteristics (nominal)

source	APS bending magnet
critical energy	19.51 keV
energy range	1.0 - 100.0 keV
on-axis peak brilliance at 16.3 keV	$2.9 \times 10^{15}$ ph/sec/mrad <sup>2</sup> /mm <sup>2</sup> /0.1%bw
on-axis peak angular flux at 16.3 keV	$9.6 \times 10^{13}$ ph/sec/mrad <sup>2</sup> /0.1%bw
on-axis peak horizontal angular flux at 5.6 keV	$1.6 \times 10^{13}$ ph/sec/mradh/0.1%bw
source size at critical energy $\Sigma_x$ $\Sigma_y$	145 $\mu\text{m}$ 36 $\mu\text{m}$
source divergence at critical energy $\Sigma_{x'}$ $\Sigma_{y'}$	6 mrad 47 $\mu\text{rad}$

## Optics & Optical Performance

- 4.5–80 keV spectral range (white beam)

## Experiment Stations

- 5-BM-B  
monochromatic beam station
- 5-BM-C  
monochromatic and white beam beam station  
powder diffraction
- 5-BM-D  
monochromatic and white beam station  
general purpose station for  
EXAFS  
crystallography  
x-ray scattering

## Detectors

- spectroscopy-grade ionization chambers
- Lytle detector
- Fuji BAS 2000 image-plate system
- scintillation detectors
- Mar 135 mm CCD detector
- Spectra Source CCD camera for imaging applications
- EG&G Ortec Iglet solid-state detector

## Beamline Controls and Data Acquisition

- Linux running custom control code for motors and actuators on all systems
- SPEC used to control diffractometers and surface-science instruments
- Vendor-supplied software used to control CCD detectors

## Beamline Support Equipment/Facilities

- 4-circle diffractometer
- 2-circle powder diffractometer
- Mar CCD detector system
- Fuji BAS 2000 image-plate system

**Beamline contacts:** John Quintana  
JerrieLea Hopf

tel 630.252.0221  
tel 630.252.0222

jpq@nwu.edu  
jlhopf@nwu.edu

# Beamline 5-ID / DND-CAT

**Scientific focus:** Crystallography, EXAFS, polymer scattering, and surface science

**Scientific programs:** Material science and engineering, polymer science, crystallography, surface science, and catalysis

## Insertion Device Source Characteristics (nominal)

source	Undulator A
period	3.30 cm
length	2.47 m
effective $K_{\max}$ (at minimum gap = 10.5 mm)	2.78
energy range 1st harmonic	2.9 - 13.0 keV
energy range 1st - 5th harmonics	2.9 - 45.0 keV
on-axis peak brilliance at 6.5 keV	$9.6 \times 10^{18}$ ph/sec/mrad <sup>2</sup> /mm <sup>2</sup> /0.1% bw
source size at 8.0 keV $\sum_x$ $\sum_y$	359 $\mu\text{m}$ 21 $\mu\text{m}$
source divergence at 8.0 keV $\sum_{x'}$ $\sum_{y'}$	24 $\mu\text{rad}$ 6.9 $\mu\text{rad}$

## Optics & Optical Performance

- 7–18 keV spectral range
- cryogenic Si(111) crystals  
fixed-gap channel cut configuration
- energy range limited due to crystal size

## Experiment Stations

- 5-ID-B  
monochromatic general-purpose beam station  
4-circle Huber goniometer  
Mar CCD  
single-crystal diffraction  
EXAFS
- 5-ID-C  
monochromatic beam station  
surface science using custom UHV system  
surface scattering  
standing waves
- 5-ID-D  
small- and wide-angle scattering  
polymer scattering station

## Detectors

- spectroscopy-grade ionization chambers
- Lytle detector
- Fuji BAS 2000 image-plate system
- scintillation detectors
- Mar 135 mm CCD detector
- Spectra Source CCD camera for imaging applications
- EG&G Ortec Iglet solid-state detector

## Beamline Controls and Data Acquisition

- Linux running custom control code for motors and actuators for all systems
- SPEC used to control diffractometers and surface-science instruments
- Vendor-supplied software used to control CCD detectors

## Beamline Support Equipment/Facilities

- 6-circle diffractometer
- 2-circle powder diffractometer
- custom UHV system
- Mar CCD detector system
- Fuji BAS 2000 image-plate system

**Beamline contacts:** John Quintana  
JerrieLea Hopf

tel 630.252.0221  
tel 630.252.0222

jpq@nwu.edu  
jlhopf@nwu.edu





## Midwest Universities CAT

---

### Overview

The Midwest Universities CAT was organized in 1990. At present, its member institutions include Ames Laboratory/Iowa State University, University of Missouri at Columbia, Georgia Institute of Technology, Washington University, University of Wisconsin at Madison, Kent State University, State University of New York at Stony Brook, Michigan State University, and the Institut für Festkörperforschung Forschungszentrum Jülich GmbH. The main undulator line of the sector features a four-circle diffractometer for single-crystal work, a liquid surface reflectometer, and a new surface science chamber. In addition, a high-energy side station (30-120 keV) is currently under construction and will run in parallel with the main undulator line [3-22 keV for Si (111) monochromating crystals].

### Research Focus

Research efforts are focused in the areas of magnetic scattering and spectroscopy, as well as in resonant- and nonresonant-scattering studies of magnetic materials. The technique of circular magnetic x-ray dichroism (CMXD) offers tremendous potential as a tool for studies of the electronic states associated with magnetic order in ferromagnetic materials of fundamental and technological interest. Resonant and nonresonant magnetic x-ray scattering measurements offer important and complementary means of determining magnetic structures in materials that are ill-suited, by reasons of size or chemical composition, to traditional neutron measurements (e.g., for magnetic structure studies of metastable phases and surfaces).

Research efforts in the surface-scattering program are centered on the study of the kinetics and growth of two-dimensional systems, the role of defects in epitaxy, ordered non-epitaxial overlayers, phase transitions, and investigations of liquid surfaces. The surface chamber constructed for this sector will enable surface-scattering investigations to expand into new areas: real-time growth characterization in a variety of growth environments, low-Z adsorbate systems, and two-dimensional surface kinetics. A liquid surface diffractometer will be used to probe the chemistry and physics of monolayer films at liquid surfaces, as well as of realistic models of biological membranes.

General scattering programs include microdiffraction techniques, coherent scattering, high-temperature processing of materials, studies of liquid crystals and polymers, and strain in thin films and interfaces.

Work on the high-energy side station will include techniques such as nonresonant magnetic scattering, radial distribution function analysis of disordered and partially ordered materials, and *in situ* time-resolved measurements of high-temperature materials processing.

---

<b>CAT contacts:</b>	Alan Goldman, Director	tel 515.294.3585	goldman@ameslab.gov
	Paul Miceli, Deputy Dir.	tel 573.882.8328	pfm@wald.physics.missouri.edu
	Doug Robinson, Beamline Scientist	tel 630.252.0247	drobinsn@iastate.edu
	Laura Morisco, CAT Coordinator	tel 630.252.0243	

# Beamline 6-ID / $\mu$ -CAT

**Scientific focus:** Materials science

**Scientific programs:** Magnetic scattering and spectroscopy, liquid surface studies, and surface scattering

## Insertion Device Source Characteristics (nominal)

source	Undulator A
period	3.30 cm
length	2.47 m
effective $K_{\max}$ (at minimum gap = 10.5 mm)	2.78
energy range 1st harmonic	2.9 - 13.0 keV
energy range 1st - 5th harmonics	2.9 - 45.0 keV
on-axis peak brilliance at 6.5 keV	$9.6 \times 10^{18}$ ph/sec/mrad <sup>2</sup> /mm <sup>2</sup> /0.1% bw
source size at 8.0 keV $\sum_x$ $\sum_y$	359 $\mu$ m 21 $\mu$ m
source divergence at 8.0 keV $\sum_{x'}$ $\sum_{y'}$	24 $\mu$ rad 6.9 $\mu$ rad

## Optics & Optical Performance

- Kohzu double-crystal constant-offset monochromator
  - 29.5 m from source
  - 3.2–20 keV energy range Si(111)
  - $10^{-4}$  energy resolution ( $\Delta E/E$ ) at 10 keV
  - 25 mm offset
  - liquid-nitrogen cooling
- 3-stripe vertically focusing 70 cm mirror
  - 60 cm active length (clear aperture)
  - Al, Rh, Pt coating stripes
  - 12, 23, 31 keV high-energy cutoffs at  $0.15^\circ$
  - 1:1 vertical focusing (variable) at 6-ID-B
  - < 0.5 arcsec meridional slope error
  - < 5 arcsec sagittal slope error
  - 3 Å rms roughness

## Experiment Stations

- 6-ID-A
  - white beam first optics enclosure
  - Kohzu monochromator
- 6-ID-B
  - monochromatic beam station
  - 4-circle diffractometer
  - liquid surface reflectometer
- 6-ID-C
  - monochromatic beam station
  - UHV surface science chamber

## Detectors

- NaI detectors
- ionization chambers
- 3-element Si energy dispersive array detector

## Beamline Controls and Data Acquisition

- Sun UNIX running EPICS with VME
- SPEC control software

## Beamline Support Equipment/Facilities

- 6K Displex
- polarization modifier
- polarization analyzer
- high-temperature furnace for powder diffraction

**Beamline contacts:** Alan Goldman  
Doug Robinson

tel 515.294.3585  
tel 630.252.0247

goldman@ameslab.gov  
drobinson@iastate.edu



## MHATT-CAT

### Center for Real-Time X-ray Studies

---

#### Overview

The Center for Real-Time X-ray Studies CAT (MHATT-CAT) was formed in 1989 as a university-industry partnership between the University of Michigan, Howard University, and Bell Laboratories (Lucent Technologies). MHATT-CAT's facilities include two fully instrumented beamlines (one based on an undulator and the other on a bend magnet) to exploit the unique characteristics of the APS. A high-speed communications network between the participating institutions and the APS is planned to facilitate the involvement of students in the MHATT-CAT collaboration and to maximize participation in research at national user facilities.

#### Research Focus

MHATT-CAT's research program is focused on time-dependent structural phenomena across a broad range of materials systems, including thin films, polymers, liquid crystals, and biomaterials. The undulator line serves as the host for microprobe studies of high-temperature superconductors, interconnects, and optical fibers; for the development of coherent x-ray scattering techniques; and for the application of ultrafast lasers. The sector is expected to develop into an international gathering point for scientists and engineers with strong interest in time-resolved structural studies.

---

**CAT contact:**

Roy Clarke, Director  
WWW Site

tel 734.764.4466

royc@umich.edu  
<http://www.mhatt.aps.anl.gov>

# Beamline 7-ID / MHATT-CAT

**Scientific focus:** Real-time x-ray studies

**Scientific programs:** Real-time x-ray studies, x-ray microcharacterization, and x-ray photon correlation spectroscopy

## Insertion Device Source Characteristics (nominal)

source	Undulator A
period	3.30 cm
length	2.47 m
effective $K_{\max}$ (at minimum gap = 10.5 mm)	2.78
energy range 1st harmonic	2.9 - 13.0 keV
energy range 1st - 5th harmonics	2.9 - 45.0 keV
on-axis peak brilliance at 6.5 keV	$9.6 \times 10^{18}$ ph/sec/mrad <sup>2</sup> /mm <sup>2</sup> /0.1% bw
source size at 8.0 keV $\sum_x$ $\sum_y$	359 $\mu\text{m}$ 21 $\mu\text{m}$
source divergence at 8.0 keV $\sum_{x'}$ $\sum_{y'}$	24 $\mu\text{rad}$ 6.9 $\mu\text{rad}$

## Optics & Optical Performance

- BESSRC-type double-crystal fixed-exit HHL monochromator
  - 30 m from source
  - 5–50 keV energy range
  - $10^{-4}$  energy resolution ( $\Delta E/E$ ) at 10 keV
  - 35 mm offset
  - liquid-nitrogen cooling
- pink-beam mirror filter
  - 30 m from source (centered on HHL mono)
  - 6–10 keV energy range
  - 0.03 energy resolution ( $\Delta E/E$ )

## Experiment Stations

- 7-ID-A
  - white beam first optics enclosure
  - HHL monochromator
- 7-ID-B
  - white beam station
  - 5 m x 7 m
- 7-ID-C
  - monochromatic and pink beam station
  - 5 m x 7 m
  - diffraction
- 7-ID-D
  - monochromatic and pink beam station
  - 5 m x 6 m
  - diffraction

## Detectors

- NaI scintillation
- ionization chambers
- various CCD detectors
- Amptek Si detector

## Beamline Controls and Data Acquisition

- EPICS
- SPEC
- k-Space Associates real-time CCD data acquisition
- LabView VI environment



## Beamline Support Equipment

- Newport/Micro-Control 6-circle Kappa goniometer (7-ID-C)
- Huber 4-circle goniometer (7-ID-B)
- two xyz optical tables (4 ft x 6 ft)

**Beamline contacts:** Eric Dufresne

tel 630.252.0274

dufresne@umich.edu

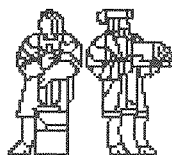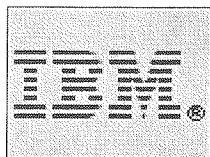

## IBM-MIT-McGill-Whitehead CAT

---

### Overview

IBM-MIT-McGill-Whitehead CAT (IMMW-CAT) was organized by individual scientists to build and operate a sector at the APS optimized to conduct many kinds of x-ray scattering research. The undulator beamline is well suited to small-angle scattering experiments; high-resolution, wide-angle scattering experiments; and experiments involving dynamic phenomena in materials science, condensed matter physics, and biological physics. The bending magnet beamline is being developed for macromolecular crystallography.

### Research Focus

IMMW-CAT's research spans a broad range of topics. Experiments conducted by CAT-member scientists fall primarily into the following areas:

1. Small-angle x-ray scattering studies, especially of soft matter.
2. Time-resolved x-ray scattering studies of bulk phase transition kinetics.
3. High-resolution, wide-angle diffraction studies of solid state materials, their phase behavior, and phase transitions.
4. Studies of equilibrium and non-equilibrium dynamics using the emerging technique of x-ray photon correlation spectroscopy (XPCS) and the development of XPCS-related coherent x-ray methods.
5. High-resolution x-ray diffraction studies of the structure of biologically important macromolecules (including MAD studies).

---

<b>CAT contacts:</b>	Simon Mochrie, Director	tel 617.253.6588	simon@lindy.mit.edu
	D. Mark Sutton, Deputy Dir.	tel 514.398.6523	mark@physics.mcgill.ca

# Beamline 8-ID / IMMW-CAT

**Scientific focus:** Condensed matter physics and materials science

**Scientific programs:** Intensity fluctuation spectroscopy studies using coherent x-rays, time-resolved x-ray scattering, x-ray scattering studies and very low temperatures, and x-ray magnetic scattering

## Insertion Device Source Characteristics (nominal)

source	Undulator A
period	3.30 cm
length	2.47 m
effective $K_{\max}$ (at minimum gap = 10.5 mm)	2.78
energy range 1st harmonic	2.9 - 13.0 keV
energy range 1st - 5th harmonics	2.9 - 45.0 keV
on-axis peak brilliance at 6.5 keV	$9.6 \times 10^{18}$ ph/sec/mrad <sup>2</sup> /mm <sup>2</sup> /0.1% bw
source size at 8.0 keV $\sum_x$ $\sum_y$	359 $\mu\text{m}$ 21 $\mu\text{m}$
source divergence at 8.0 keV $\sum_{x'}$ $\sum_{y'}$	24 $\mu\text{rad}$ 6.9 $\mu\text{rad}$

## Optics & Optical Performance

- horizontally deflecting mirror  
29.2 m from source  
6–30 keV energy range  
water cooling
- transmission diamond monochromator  
51.6 m from source  
6–30 keV energy range  
 $\Delta\lambda/\lambda=44 \times 10^{-6}$  energy resolution at 8 keV  
symmetric Bragg diamond (111)  
water cooling

## Experiment Stations

- 8-ID-A  
white beam first optics enclosure
- 8-ID-D  
pink beam monochromator enclosure
- 8-ID-E  
pink beam experimental station
- 8-ID-I  
pink beam experimental end station

## Detectors

- two PI direct detection CCD detectors
- PI optically coupled CCD detector
- photodiode array linear detector
- Amptek CZT and Si detectors

## Beamline Controls and Data Acquisition

- Linux-based control system
- SPEC x-ray control software

## Beamline Support Equipment/Facilities

- coherent small-angle scattering spectrometer
- time-resolved diffractometer
- He-3 cryostat on horizontally scattering spectrometer

## Beamline contacts:

Larry Lurio  
Simon Mochrie  
Mark Sutton

tel 630.252.0281  
tel 617.253.6588  
tel 514.398.6523

larry@slate.imm.aps.anl.gov  
simon@lindy.mit.edu  
mark@physics.mcgill.ca

## CMC-CAT

### Complex Materials Consortium CAT

---

#### Overview

The Complex Materials Consortium CAT (CMC-CAT) was organized by individuals who share a common research interest in complex materials and a strong dedication to the use of the unique properties of synchrotron radiation to structurally and functionally characterize them. This organization was subsequently formalized into a consortium containing the following institutional members: Exxon Research and Engineering Company, University of Pennsylvania, Brookhaven National Laboratory, Princeton University, University of California at Santa Barbara, Los Alamos National Laboratory, Oak Ridge National Laboratory, and the University of Tennessee.

#### Research Focus

Complex materials, broadly defined, are materials that possess unique and novel properties by virtue of the complexity of their molecules and/or the complexity of interaction between their molecules. These materials include polymers, surfactants, liquid crystals, biomaterials, membranes, and thin molecular films of hydrocarbons on solid or liquid surfaces. They also include novel synthetic materials such as fullerenes, fibers, polymer composites, oxides, and zeolites. A basic understanding of the molecular structure, morphology, and molecular dynamics of such materials, particularly under non-equilibrium conditions or at interfaces, is essential for the development and optimization of novel materials and processes for industry. It is equally important for obtaining a complete theoretical description of the equilibrium and non-equilibrium behavior of complex, multicomponent systems.

At the APS, CMC-CAT uses several experimental techniques to study complex materials on a variety of length and time scales. The techniques include small-angle and wide-angle x-ray scattering, liquid surface spectrometry, x-ray surface scattering, x-ray spectroscopy, magnetic scattering, inelastic scattering, and x-ray microtomography and imaging. The studies depend on novel technical developments such as ultrafast-high-resolution two-dimensional imaging detectors, high-speed data inversion and image-processing algorithms, and novel x-ray microfocusing devices. Using such techniques, CMC-CAT researchers study phenomena such as the dynamic response of polymers, membranes, fibers, and fluids to applied stresses; the kinetics of phase separation in mixed systems; the formation of pitting, corrosion, and protective layers on metals in electrolytic solutions; the formation of novel surface structures on catalysts undergoing reactions; the formation of cracks and voids in composite materials; the properties of transition metal oxides; and the flow of fluids through microporous solids such as oil-bearing rocks or sand.

---

<b>CAT contacts:</b>	Doon Gibbs, Director	tel 631.344.4608	doon@solids.phy.bnl.gov
	Thomas Gog, Assoc. Dir.	tel 630.252.0320	gog@anl.gov
	Chitra Venkataraman, Beamline Scientist	tel 630.252.0327	chitra@anl.gov

# Beamline 9-BM / CMC-CAT

**Scientific focus:** Materials science

**Scientific programs:** EXAFS of complex materials, x-ray surface scattering, and general diffraction

## Bending Magnet Source Characteristics (nominal)

source	APS bending magnet
critical energy	19.51 keV
energy range	2.0 - 40.0 keV
on-axis peak brilliance at 16.3 keV	$2.9 \times 10^{15}$ ph/sec/mrad <sup>2</sup> /mm <sup>2</sup> /0.1%bw
on-axis peak angular flux at 16.3 keV	$9.6 \times 10^{13}$ ph/sec/mrad <sup>2</sup> /0.1%bw
on-axis peak horizontal angular flux at 5.6 keV	$1.6 \times 10^{13}$ ph/sec/mradh/0.1%bw
source size at critical energy $\sum_x$ $\sum_y$	145 $\mu$ m 36 $\mu$ m
source divergence at critical energy $\sum_{x'}$ $\sum_{y'}$	6 mrad 47 $\mu$ rad

## Optics & Optical Performance

- white beam slits (APS design L3-20)
- double-crystal, water-cooled Si monochromator
- toroidal mirror
- monochromatic slits (APS design L4)

## Experiment Stations

- 9-BM-A  
white beam first optics enclosure
- 9-BM-B  
monochromatic beam station  
EXAFS
- 9-BM-C  
monochromatic beam station  
surface diffraction  
general diffraction

## Detectors

- Smart 1500 CCD
- scintillation counters

## Beamline Controls and Data Acquisition

- Sun UNIX running EPICS with VME, SPEC

## Beamline Support Equipment/Facilities

- UHV surface apparatus
- EXAFS instrumentation

**Beamline contacts:** Thomas Gog  
Chitra Venkataraman

tel 630.252.0320  
tel 630.252.0327

gog@anl.gov  
chitra@anl.gov

# Beamline 9-ID / CMC-CAT

**Scientific focus:** Materials science

**Scientific programs:** SAXS from complex materials, liquid/solid surface scattering, EXAFS, general diffraction, and magnetic scattering

## Insertion Device Source Characteristics (nominal)

source	Undulator A
period	3.30 cm
length	2.47 m
effective $K_{\text{max}}$ (at minimum gap = 10.5 mm)	2.78
energy range 1st harmonic	4.0 - 22.0 keV
energy range 1st - 5th harmonics	4.0 - 45.0 keV
on-axis peak brilliance at 6.5 keV	$9.6 \times 10^{18}$ ph/sec/mrad <sup>2</sup> /mm <sup>2</sup> /0.1% bw
source size at 8.0 keV $\sum_x$ $\sum_y$	359 $\mu\text{m}$ 21 $\mu\text{m}$
source divergence at 8.0 keV $\sum_{x'}$ $\sum_{y'}$	24 $\mu\text{rad}$ 6.9 $\mu\text{rad}$

## Optics & Optical Performance

Primary branch optics and optical performance

- white-beam slit
  - 4.5 mm x 4.5 mm input aperture
  - 2.1 mm x 2.1 mm max. output aperture
  - ctr. of output aperture adjustable w/ input aperture
- double-crystal Si(111) monochromator
  - 30.5 m from source
  - 3.1–22.5 keV energy range
  - 25.0 mm beam offset
  - $10^{-4}$  energy resolution ( $\Delta E/E$ ) at 10 keV
  - cryogenic cooling
- vertical-focusing mirror
  - 34.4 m from source
  - flat figure with bender
  - 600 mm x 120 mm
  - Au, Ni, Rh coatings
  - 0–4 mrad incident angle

- horizontal-focusing mirror
  - 33.2 m from source
  - flat figure with bender
  - 900 mm x 120 mm
  - Au, Ni, Rh coatings
  - 0–4 mrad incident angle

Secondary branch optics and optical performance

- side-station monochromator
  - 28 m from source
  - horizontal double crystal
  - 5–9.5 keV energy range
  - water-cooled diamond crystal
  - 1000 mm horizontal beam offset
  - $10^{-4}$  energy resolution
- vertical-focusing mirror
- horizontal-focusing mirror

## Experiment Stations

- 9-ID-A
  - white beam first optics enclosure
- 9-ID-B & -C
  - monochromatic beam station
  - general diffraction
  - SAXS
  - EXAFS
  - liquids

## Detectors

- Smart 1500 CCD
- BNL 2D WIRE
- BNL 1D WIRE

## Beamline Controls and Data Acquisition

- Sun UNIX running EPICS w/ VME, SPEC, etc.

## Beamline Support Equipment/Facilities

- 6-circle diffractometer
- SAXS set up
- Bonse-Hart camera
- liquid surface reflectometer
- UHV surface apparatus


**Beamline contacts:** Thomas Gog  
Chitra Venkataraman

tel 630.252.0320  
tel 630.252.0327

gog@anl.gov  
chitra@anl.gov





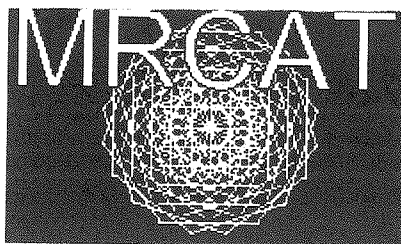

## Materials Research CAT

---

### Overview

The Materials Research CAT (MR-CAT) is composed of groups from four universities (University of Florida, Illinois Institute of Technology, Northwestern University, and University of Notre Dame), a major corporation (Beyond Petroleum) and a federal research laboratory (Argonne National Laboratory/ Chemical Technology and Environmental Research Divisions).

### Research Focus

The scientific program of MR-CAT focuses on the study of advanced materials *in situ* as a means of characterizing their structure and electronic properties, as well as understanding their preparation. The primary research techniques of wide- and small-angle scattering, diffraction (single-crystal and powder), absorption spectroscopy (XFS), reflectivity, standing waves, diffraction anomalous fine structure (DAFS), and time-dependent and microfocus techniques are being used to study the following *in situ*:

- structural phase changes, especially in non-equilibrium systems
- disordered systems (e.g., alloys and amorphous materials)
- growth, recrystallization, surfaces, and interfaces of electronic materials
- catalysts (*in situ* time-dependent studies)
- structure of static and dynamic confined liquids
- organic thin films and self-assembling systems
- polymers (e.g., dynamics of block copolymers, single-fiber studies)

---

<b>CAT contacts:</b>	Bruce Bunker, Director	tel 219.631.7219	bunker.1@nd.edu
	Carlo Segre, Deputy Dir.	tel 312.567.3498	segre@iit.edu

# Beamline 10-ID / MR-CAT

**Scientific focus:** Materials science and environmental science

**Scientific programs:** Spectroscopy of dilute environmental systems, combined techniques for *in situ* characterization of materials, and zone plate micro-focusing spectroscopy and imaging

## Insertion Device Source Characteristics (nominal)

source	Undulator A
period	3.30 cm
length	2.47 m
effective $K_{\max}$ (at minimum gap = 10.5 mm)	2.78
energy range 1st harmonic	2.9 - 13.0 keV
energy range 1st - 5th harmonics	2.9 - 45.0 keV
on-axis peak brilliance at 6.5 keV	$9.6 \times 10^{18}$ ph/sec/mrad <sup>2</sup> /mm <sup>2</sup> /0.1% bw
source size at 8.0 keV $\Sigma_x$ $\Sigma_y$	359 $\mu$ m 21 $\mu$ m
source divergence at 8.0 keV $\Sigma_{x'}$ $\Sigma_{y'}$	24 $\mu$ rad 6.9 $\mu$ rad

## Optics & Optical Performance

- Daresbury/Fisons monochromator
  - 30 m from source
  - 4–30 keV energy range Si(111)
  - 1st crystal, IIT-design, cryo cooling w/  
pseudo-fixed offset crystal cage
  - 2nd crystal, 240 mm with Bragg normal  
motion
  - $10^{-4}$  energy resolution ( $\Delta E/E$ ) at 10 keV
  - 35 mm offset
  - liquid-nitrogen cooling
- harmonic rejection mirror
- zone plate focusing

## Experiment Stations

- 10-ID-A
  - white beam first optics enclosure
  - monochromator
- 10-ID-B
  - monochromatic beam station
  - all experiments 4 m x 12 m

## Detectors

- Bruker 1k x 1k CCD
- Canberra Si energy-dispersive detector
- Daresbury/Lesker spectroscopy ionization  
chambers
- Lytle detectors
- 13-element Canberra/XIA Ge detector system

## Beamline Controls and Data Acquisition

- Debian Linux running “MX”

## Beamline Support Equipment/Facilities

- large Huber 8-circle goniometer with positioning
- X95 rail system for experiment alignment

**Beamline contacts:** Bruce Bunker  
Carlo Segre

tel 219.631.7219  
tel 312-567-3498

bunker.1@nd.edu  
segre@iit.edu



## Basic Energy Sciences Synchrotron Radiation Center CAT

---

### Overview

The Basic Energy Sciences Synchrotron Radiation Center CAT (BESSRC-CAT) includes scientists from the Chemistry and Materials Science Divisions at Argonne National Laboratory (ANL), researchers from the ANL geophysics program, and scientists from the Department of Physics at Northern Illinois University. Its two sectors are developed and instrumented to serve the research needs of users with particular interests in materials science, chemical science, atomic physics, solid-state physics, and geosciences. The facility's two sectors are equipped with seven separate instrumental stations. One sector consists of an undulator and bending magnet with four experiment stations. The other sector has an elliptical multipole wiggler and is split into three beamlines with three experiment stations.

### Research Focus

Scientific areas of interest include time-resolved studies of photoexcited states in photosynthetic materials, real-time investigations of chemical reactions, and time-dependent studies of phase transformations in crystalline compounds. Structural studies of actinide materials, studies of ultra-small crystals, trace-element analysis, and surface and interface studies also form part of the research focus. Examples of studies conducted include radiation-induced surface modification, impurity clustering at grain boundaries, measurements of the electrical double layer at electrodes in electrochemical cells, concentration profiling at dissimilar interfaces, studies of the magnetic structure of ultrathin films, and others.

BESSRC-CAT has an MOCVD system in the undulator beamline for *in situ* growth studies of oxides and GaN and also an MBE system on line for structural studies of surfaces and interfaces under UHV conditions. There is an extensive research program in the use of small-angle x-ray scattering and anomalous small-angle x-ray scattering with particular emphasis on time-resolved measurements. BESSRC-CAT has established a research program directed towards the use of high-energy diffraction techniques to study the structure of materials. Additionally, an extensive research program is directed towards the use of circularly polarized x-rays in the study of the magnetic properties of materials, i.e., using XMCD, magnetic Compton scattering, and magnetic scattering.

---

<b>CAT contacts:</b>	Pedro Montano, Director	tel 630.252.1553	montano@anl.gov
	Mark Beno, Deputy Dir.	tel 630.252.3507	beno@anl.gov
	Jeanne Cowan, Admin. Asst.	tel 630.252.1553	cowan@anl.gov

# Beamline 11-ID / BESSRC-CAT

**Scientific focus:** Spectroscopy, scattering, high-energy diffraction, and magnetic Compton scattering

**Scientific programs:** Materials science (x-ray scattering, high-energy diffraction, dichroism, magnetism), geoscience, atomic physics, and chemistry (spectroscopy [EXAFS, XANES], time-dependent spectroscopy of excited states)

## Insertion Device Source Characteristics (nominal)

source	elliptical multipole wiggler
period	16.0 cm
length	2.8 m
number of poles	36, electromagnets 34, permanent magnets
effective $K_x$ max (at max. current = 1.15 kA, $B_{max} = 0.087$ Tesla)	1.3
peak $K_y$ max (at minimum gap = 24.0 mm)	14.4
switching frequency	0 - 10 Hz
critical energy (at minimum gap = 24.0 mm)	59.1 keV
energy range	5.0 - 200.0 keV
source size at critical energy $\sum_x$ $\sum_y$	359 $\mu\text{m}$ 21 $\mu\text{m}$
source divergence at critical energy $\sum_{x'}$ (FWHM 1.9 mrad, non-Gaussian) $\sum_{y'}$	820 $\mu\text{rad}$ 47 $\mu\text{rad}$

## Optics & Optical Performance

- focusing Bragg or flat Laue monochromator
- 11-ID-D
  - modified standard BESSRC monochromator
  - 4.8–50 keV energy range Si(220)
  - 20 mm fixed offset
  - cryo-cooled 1st and 2nd crystals
  - mirrors
    - Pt and Pd coatings
    - 2.8 mrad incident angle/Pt mirror
    - 4.0 mrad incident angle/Pd mirror
- 11-ID-B & -C
  - Bragg or Laue horiz. deflecting monochromator
  - 40–200 keV energy range
  - Si(311) or Si(220) crystal
  - water/Ga cooling

## Experiment Stations

- 11-ID-A
  - white beam first optics enclosure
  - double-crystal monochromator (for 11-ID-D)
  - double (side-by-side) mirror system (for 11-ID-D)
  - side-deflecting monochromator (for B & C stations)
- 11-ID-B
  - magnetic Compton scattering
- 11-ID-C
  - high-energy diffractometer
- 11-ID-D
  - spectroscopy (EXAFS, XANES, dichroism)
  - x-ray scattering

## Detectors

- ionization chambers
- Bicron
- solid-state detectors (Si and Ge)
- Canberra 9-element Ge

## Beamline Controls and Data Acquisition

- 11-ID-B & -C
  - MacOS-EPICS
- 11-ID-B
  - Solaris-EPICS
  - SPEC
  - MacOS
  - Windows EPICS

## Beamline Support Equipment/Facilities

- 11-ID-C
  - high-energy diffractometer (vertical, post mono)
  - diffractometer
  - analyzer
  - detector
- 11-ID-D
  - spectroscopy table
  - 8-circle Huber  $\psi$  goniometer

**Beamline contacts:** Anneli Munkholm  
Uta Ruett  
Mark Beno

tel 630.252.0395  
tel 630.252.0383  
tel 630.252.3507

munkholm@anl.gov  
ruett@msd.anl.gov  
beno@anl.gov

# Beamline 12-BM / BESSRC-CAT

**Scientific focus:** Spectroscopy, x-ray scattering, and diffraction

**Scientific programs:** Materials science (x-ray scattering, diffraction, spectroscopy), geoscience (mineral-liquid interface characterization by x-ray scattering and spectroscopy), chemistry (EXAFS and XANES spectroscopy of actinides and lanthanides), and atomic physics

## Bending Magnet Source Characteristics (nominal)

source	APS bending magnet
critical energy	19.51 keV
energy range	1.0 - 100.0 keV
on-axis peak brilliance at 16.3 keV	$2.9 \times 10^{15}$ ph/sec/mrad <sup>2</sup> /mm <sup>2</sup> /0.1%bw
on-axis peak angular flux at 16.3 keV	$9.6 \times 10^{13}$ ph/sec/mrad <sup>2</sup> /0.1%bw
on-axis peak horizontal angular flux at 5.6 keV	$1.6 \times 10^{13}$ ph/sec/mradh/0.1%bw
source size at critical energy $\sum_x$ $\sum_y$	145 $\mu$ m 36 $\mu$ m
source divergence at critical energy $\sum_{x'}$ $\sum_{y'}$	6 mrad 47 $\mu$ rad

## Optics & Optical Performance

- BESSRC standard monochromator  
4–34 keV energy range Si(111)  
35 mm offset  
water cooling  
UHV vacuum compatibility
- Mirrors (removable double-mirror system)  
down deflection, flat figure, Pd coating  
up deflection, toroidal figure, Rh coating  
20 keV energy cutoff

## Experiment Stations

- 12-BM-A  
white beam first optics enclosure  
white beam slits  
monochromator and mirror system
- 12-BM-B  
monochromatic beam station  
x-ray scattering  
diffraction  
spectroscopy (EXAFS and XANES)

## Detectors

- ionization chambers
- Bicron
- Lytle detector
- solid-state detectors
- single-element Si and Ge
- Canberra 9-element Ge

## Beamline Controls and Data Acquisition

- SPARC station with Solaris operating system running EPICS with VME and SPEC software
- Windows NT running EPICS applications
- Macintosh with MacOS software running EPICS applications

## Beamline Support Equipment/Facilities

- 12-BM-B  
spectroscopy table  
6-circle Huber goniometer

**Beamline contacts:** Jennifer Linton  
Mark Beno  
Pedro Montano

tel 630.252.5007  
tel 630.252.3507  
tel 630.252.6239

linton@msd.anl  
beno@anl.gov  
montano@anl.gov

# Beamline 12-ID / BESSRC-CAT

**Scientific focus:** Materials science, atomic physics, and geoscience

**Scientific programs:** Materials science (elastic and inelastic scattering, MOCVD, x-ray standing waves, electrochemistry), geoscience (mineral–liquid interface characterization), chemistry (small-angle scattering), and atomic physics (spectroscopy)

## Insertion Device Source Characteristics (nominal)

source	Undulator A
period	3.30 cm
length	2.47 m
effective $K_{\max}$ (at minimum gap = 10.5 mm)	2.78
energy range 1st harmonic	2.9 - 13.0 keV
energy range 1st - 5th harmonics	2.9 - 45.0 keV
on-axis peak brilliance at 6.5 keV	$9.6 \times 10^{18}$ ph/sec/mrad <sup>2</sup> /mm <sup>2</sup> /0.1% bw
source size at 8.0 keV $\sum_x$ $\sum_y$	359 $\mu\text{m}$ 21 $\mu\text{m}$
source divergence at 8.0 keV $\sum_{x'}$ $\sum_{y'}$	24 $\mu\text{rad}$ 6.9 $\mu\text{rad}$

## Optics & Optical Performance

- BESSRC standard monochromator
  - Si(111) crystal
  - 4–27 keV energy range
  - 35 mm offset
  - liquid-nitrogen cooling
  - 2nd crystal: Cu braid, liquid nitrogen cooling
  - UHV vacuum compatibility
- flat focusing monochromatic mirror
  - 40 m from source
  - Pt, Pd, and SiO<sub>2</sub> coatings
  - 2.5–4 mrad incident angle
  - ~100  $\mu\text{m}$  focus in 12-ID-B

## Experiment Stations

- 12-ID-A
  - white beam first optics enclosure
  - monochromator
  - quad photodiode
  - monochromator pinhole apertures

- 12-ID-B
  - spectroscopy
  - elastic and inelastic scattering
  - reflectivity
  - small-angle scattering
- 12-ID-C
  - atomic physics
  - small-angle scattering (low Q)
- 12-ID-D
  - MOCVD chamber
  - surface chamber
  - x-ray scattering

## Detectors

- ionization chambers
- Bicron
- solid-state detectors
- single-element Si and Ge
- Canberra 9-element Ge
- 2-D proportional wire detector
- gold CCD camera

## Beamline Controls and Data Acquisition

- 12-ID-B
  - Solaris running SPEC, Windows NT, and MacOS running EPICS applications
- 12-ID-C
  - Solaris running SPEC, Windows NT, and MacOS running EPICS applications
- 12-ID-D
  - Linux workstation running SPEC (non-EPICS version) and EPICS applications


## Beamline Support Equipment/Facilities

- 12-ID-B
  - Huber 8-circle goniometer
  - spectroscopy
  - optics table
  - small-angle scattering instrumentation
- 12-ID-D
  - 6-circle goniometer
  - surface chamber
  - MOCVD apparatus

**Beamline contacts:** Jennifer Linton  
Mark Beno  
Pedro Montano

tel 630.252.5007  
tel 630.252.3507  
tel 630.252.6239

linton@msd.anl.gov  
beno@anl.gov  
montano@anl.gov

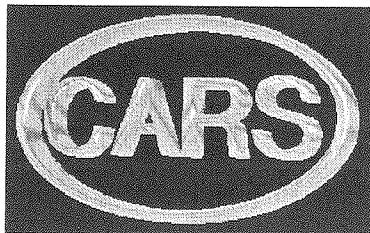

## **Consortium for Advanced Radiation Sources CAT**

---

### **Overview**

The Consortium for Advanced Radiation Sources (CARS) currently has four institutional members—The University of Chicago (the managing agent for CARS), Northern Illinois University, Southern Illinois University, and the Australian Nuclear Science and Technology Organization—and four national user groups—BioCARS for structural biology, GeoCARS for geophysical sciences, SoilEnviroCARS for soil/environmental sciences, and ChemMatCARS for chemistry and materials science. At present, CARS has developed or is developing five x-ray beamlines for use by the faculty and staff of its institutional members and national user groups for frontier research that exploits the unique characteristics of the APS as an x-ray source.

### **Research Focus**

The scientific program and experiments carried out on CARS beamlines exploit the unique features of the x-rays emitted by both wiggler and undulator sources. Two main themes are 1) the study of samples with which x-rays interact very weakly because of their small size (area or volume), low concentration, large crystallographic unit cell dimensions, or necessity for monochromatic radiation of unusually narrow bandpass and 2) time-dependent studies of samples with which x-rays interact more strongly.

For BioCARS, the main thrust is to understand basic biological processes in structural terms. Experiments involve crystallographic studies of viruses, ribosomes, and other complexes with very large unit cells; studies of microcrystals; time-resolved crystallography; and phase determination by MAD techniques.

In GeoSoilEnviroCARS, research is designed to further knowledge of the composition, structure, and properties of earth and planetary materials and the processes they control. High-pressure research in both a diamond-anvil cell and a large-volume press; x-ray diffraction, scattering, and absorption spectroscopies from earth and planetary materials; and x-ray fluorescence microprobe analysis and microtomography should make possible major advances in understanding such diverse processes as earthquake development and nitrogen fixation in roots.

In ChemMatCARS, research focuses on several aspects of dynamic and structural condensed matter and materials chemistry. Research areas include surface and interfacial properties in soft condensed matter; molecular liquids and liquid metals; chemical crystallography with emphasis on dynamics of energy transfer, fast photochemical processes, and valence and element specific crystallography; the structure and properties of molecular aggregates; and interfacial and bulk properties of novel polymers and composites including supramolecular, mesoscopic, and nanometer structures.

---

<b>CAT contacts:</b>	P. James Viccaro, Acting CAT Dir.	tel 630.252.0464	viccaro@cars.uchicago.edu
BioCARS	Keith Moffat, Assoc. Dir.	tel 773.702.9950	moffat@cars.uchicago.edu
	Wilfried Schildkamp, Assoc. Dir.	tel 630.252.0445	wilfried@cars.uchicago.edu
ChemMatCARS	P. James Viccaro, Assoc. Dir.	tel 630.252.0464	viccaro@cars.uchicago.edu
GSECARS	Mark Rivers, Assoc. Dir.	tel 630.252.0422	rivers@cars.uchicago.edu
	Steve Sutton, Assoc. Dir.	tel 630.252.0426	sutton@cars.uchicago.edu
	Joy Talsma, Exec. Admin.	tel 773.702.9506	talsma@cars.uchicago.edu

# Beamline 13-BM / GSECARS-CAT

**Scientific focus:** Geosciences, environmental science, and soil science

**Scientific programs:** High-pressure diffraction in diamond-anvil cell, high-pressure diffraction in multi-anvil press, microspectroscopy and XRF microprobe, and microtomography

## Bending Magnet Source Characteristics (nominal)

source	APS bending magnet
critical energy	19.51 keV
energy range	1.0 - 100.0 keV
on-axis peak brilliance at 16.3 keV	$2.9 \times 10^{15}$ ph/sec/mrad <sup>2</sup> /mm <sup>2</sup> /0.1%bw
on-axis peak angular flux at 16.3 keV	$9.6 \times 10^{13}$ ph/sec/mrad/0.1%bw
on-axis peak horizontal angular flux at 5.6 keV	$1.6 \times 10^{13}$ ph/sec/mradh/0.1%bw
source size at critical energy $\sum_x$ $\sum_y$	145 $\mu$ m 36 $\mu$ m
source divergence at critical energy $\sum_{x'}$ $\sum_{y'}$	6 mrad 47 $\mu$ rad

## Optics & Optical Performance

- VG double-crystal monochromator
  - 28 m from source
  - Si(111) or Si(220) crystal
  - water cooling
  - +30 mm offset
- vertical focusing mirror
  - 45 m from source
  - Si substrate
  - Pt, Rh, Si coatings
  - internal water cooling
  - elliptically bent flat geometry
  - 0–5 mrad grazing angle
  - 6:1 to 3:1 demagnification
- small Kirkpatrick-Baez microfocusing mirrors

## Experiment Stations

- 13-BM-A
  - white beam first optics enclosure
- 13-BM-B
  - white beam second optics enclosure
- 13-BM-C
  - monochromatic beam station
  - diffraction
  - tomography
- 13-BM-D
  - white or monochromatic beam station
  - multi-anvil press high-pressure diffraction
  - tomography
  - microprobe
  - microspectroscopy

## Detectors

- Canberra 16-element Ge detectors (two)
- Bruker 2K CCD detector
- Canberra single-element Ge and Si(Li) detectors
- Princeton Instruments visible light CCD cameras (two)

## Beamline Controls and Data Acquisition

- Windows NT workstations running EPICS with VME
- SPEC, IDL, Bruker SMART and GADDS
- Princeton Instruments WinView and WinSpec

## Beamline Support Equipment/Facilities

- 250-ton multi-anvil press with DIA and T-cup tooling (13-BM-D)
- laser Raman system in support laboratory

**Beamline contacts:** Mark Rivers  
Stephen Sutton  
Nancy Lazarz

tel 630.252.0422  
tel 630.252.0426  
tel 630.252.0423

rivers@cars.uchicago.edu  
sutton@cars.uchicago.edu  
lazarz@cars.uchicago.edu



# Beamline 13-ID / GSECARS-CAT

**Scientific focus:** Geosciences, environmental science, and soil science

**Scientific programs:** High-pressure diffraction in diamond-anvil cell, high-pressure diffraction in multi-anvil press, microspectroscopy and XRF microprobe, microcrystal diffraction, surface diffraction, and x-ray absorption fine structure spectroscopy

## Insertion Device Source Characteristics (nominal)

source	Undulator A
period	3.30 cm
length	2.47 m
effective $K_{\max}$ (at minimum gap = 10.5 mm)	2.78
energy range 1st harmonic	2.9 - 13.0 keV
energy range 1st - 5th harmonics	2.9 - 45.0 keV
on-axis peak brilliance at 6.5 keV	$9.6 \times 10^{18}$ ph/sec/mrad <sup>2</sup> /mm <sup>2</sup> /0.1% bw
source size at 8.0 keV $\sum_x$ $\sum_y$	359 $\mu\text{m}$ 21 $\mu\text{m}$
source divergence at 8.0 keV $\sum_{x'}$ $\sum_{y'}$	24 $\mu\text{rad}$ 6.9 $\mu\text{rad}$

## Optics & Optical Performance

- VG double-crystal monochromator
  - 29 m from source
  - Si(220) crystal
  - cryogenic cooling
  - +50 mm offset
- vertical focusing mirror
  - 45 m from source
  - Si substrate
  - Pt, Rh, Si coatings
  - internal water cooling
  - Kirkpatrick-Baez geometry
  - 0–5 mrad grazing angle
  - 6:1 to 3:1 demagnification
- horizontal focusing mirror
  - 47 m from source
  - Si substrate
  - Pt, Rh, Si coatings
  - internal water cooling
  - Kirkpatrick-Baez geometry
  - 0–5 mrad grazing angle
  - 7:1 to 3:1 demagnification
- small Kirkpatrick-Baez microfocusing mirrors

## Experiment Stations

- 13-ID-A
  - white beam first optics enclosure
- 13-ID-B
  - white beam second optics enclosure
- 13-ID-C
  - white or monochromatic beam station
  - x-ray absorption fine-structure spectroscopy
  - XRF microprobe
  - microspectroscopy
  - microcrystal diffraction
  - surface diffraction
- 13-ID-D
  - white or monochromatic beam station
  - multi-anvil press high-pressure diffraction
  - diamond anvil cell high-pressure diffraction

## Detectors

- Canberra 16-element Ge detectors (two)
- Bruker 2K CCD detector
- Canberra single-element Ge and Si(Li) detectors
- Princeton Instrument visible light CCD cameras (two)

## Beamline Controls and Data Acquisition

- Windows NT workstations running EPICS with VME
- SPEC, IDL, Bruker SMART and GADDS
- Princeton Instruments WinView and WinSpec

## Beamline Support Equipment/Facilities

- 13-ID-C
  - Newport general-purpose diffractometer
- 13-ID-D
  - diamond-anvil cell with laser heating
  - two-circle diffractometer
  - 1000-ton multi-anvil press w/ DIA and T-cup tooling
- laser Raman system in support laboratory

**Beamline contacts:** Mark Rivers  
Stephen Sutton  
Nancy Lazarz

tel 630.252.0422  
tel 630.252.0426  
tel 630.252.0423

rivers@cars.uchicago.edu  
sutton@cars.uchicago.edu  
lazarz@cars.uchicago.edu

# Beamline 14-BM / BioCARS-CAT

**Scientific focus:** Structural biology

**Scientific programs:** Large unit cell (virus) crystallography, small unit cell (protein) crystallography, MAD phasing, time-resolved crystallography, Laue diffraction, and study of microcrystals

## Bending Magnet Source Characteristics (nominal)

source	APS bending magnet
critical energy	19.51 keV
energy range	1.0 - 100.0 keV
on-axis peak brilliance at 16.3 keV	$2.9 \times 10^{15}$ ph/sec/mrad <sup>2</sup> /mm <sup>2</sup> /0.1%bw
on-axis peak angular flux at 16.3 keV	$9.6 \times 10^{13}$ ph/sec/mrad <sup>2</sup> /0.1%bw
on-axis peak horizontal angular flux at 5.6 keV	$1.6 \times 10^8$ ph/sec/mradh/0.1%bw
source size at critical energy $\Sigma_x$ $\Sigma_y$	145 $\mu$ m 36 $\mu$ m
source divergence at critical energy $\Sigma_x'$ $\Sigma_y'$	6 mrad 47 $\mu$ rad

## Optics & Optical Performance

14-BM-C optics and optical performance

- 7–17 keV energy range
- 0.21 mm hor. x 0.35 mm vert. focal spot size
- CARS-design conical focusing mirror
  - 25.731 m from source (on orbit along 3 mrad line)
  - Si substrate, Rh coating
  - water cooling at midplane
  - 4.1 mrad design angle
  - horizontal focus: focusing monochromator
  - vertical focus: bender
- CARS-design bent Ge(111) monochromator
  - 46.6 m from source (along 3 mrad centerline)
  - Ge(111) single crystal
  - water cooling
  - bent triangle horizontal focusing
- Misc. slits, collimators, filters, etc.

14-BM-D optics and optical performance

- 6.5–18.5 keV energy range
- 0.4 mm hor. x 0.6 mm vert. focal spot size
- CARS-design Si(111) double-bounce monochromator
  - 23.860 m from source (on orbit along 5 mrad line)
  - 6.5–18.5 energy range Si(111) crystal sets
  - $10^{-4}$  energy resolution ( $\Delta E/E$ ) at 10 keV
  - 1.5 in. nominal offset (fixed-exit, up bounce)
  - water cooling
- CARS-design cylindrical focusing mirror
  - Si substrate, Rh coating
  - water cooling at midplane
  - 25.731 m from source (variable height along 5 mrad line)
  - 4.1 mrad design angle
  - horizontal focus: sagittal cylindrical figure
  - vertical focus: bender mechanism
- misc. slits, collimators, filters, etc.

## Experiment Stations

- 14-BM-A
  - white beam first optics enclosure
  - 8.6 m x 1.8 m x 2.8 m (L x W x H)
- 14-BM-B
  - pink beam optics enclosure
  - 7.6 m x 1.2 m x 2.8 m (L x W x H)
- 14-BM-C
  - monochromatic beam station
  - 5.2 m x 3.7 m x 2.8 m (L x W x H)
  - virus crystallography
- 14 BM-D
  - pink or monochromatic beam station
  - 6.4 m x 2.2 m x 2.8 m (L x W x H)
  - MAD phasing
  - time-resolved Laue
  - protein crystallography

## Detectors

- 60° kappa diffractometer (14-BM-D)
- single-axis diffractometer (14-BM-C)
- ADSC Q1, ADSC Q4, MAR345 and off-line image plate detectors
- NaI scintillation detectors
- Ge detector

*continued on pg. 41*

**Beamline contacts:** Keith Moffat  
Wilfried Schildkamp  
Gary Navrotsky

tel 773.702.9950      moffat@cars.uchicago.edu  
tel 630.252.0445      schildkamp@cars.uchicago.edu  
tel 630.252.0444      navrotsky@cars.uchicago.edu

# Beamline 14-BM / BioCARS-CAT

**Scientific focus:** Structural biology

**Scientific programs:** Large unit cell (virus) crystallography, small unit cell (protein) crystallography, MAD phasing, time-resolved crystallography, Laue diffraction, and study of microcrystals

*continued from pg. 40*

## Beamline Controls and Data Acquisition

- beamline controls: HP9000/778 HPUX 10.2 running EPICS via VME
- experiment: SGI 02 IRIX 6.3 running ADSC or MAR control software via in-house (kbscan>) software
- analysis: fast SGI 02 IRIX 6.5 running DENZO, MOSFILM, DPS, etc.

## Beamline Support Equipment/Facilities

- cryo-coolers: Oxford CryoStream, Oxford CryoJet, MSC, and CARS LN<sub>2</sub>/LHe<sub>2</sub> Cooler
- collimators, filters, slits, beam stop, CCD alignment camera
- beam position monitors (1  $\mu$ m resolution)
- beam flux monitors
- BL3 facility, sample prep areas, cold room

**Beamline contacts:** Keith Moffat  
Wilfried Schildkamp  
Gary Navrotski

tel 773.702.9950  
tel 630.252.0445  
tel 630.252.0444

moffat@cars.uchicago.edu  
schildkamp@cars.uchicago.edu  
navrotski@cars.uchicago.edu

# Beamline 14-ID / BioCARS-CAT

**Scientific focus:** Structural biology

**Scientific programs:** Large unit cell (virus) crystallography, small unit cell (protein) crystallography, MAD phasing, time-resolved crystallography, Laue diffraction, and study of microcrystals

## Insertion Device Source Characteristics (nominal)

source	wiggler A
period	8.50 cm
length	2.4 m
peak $K_{\max}$ (at minimum gap = 18.1 mm)	8.74
critical energy (at minimum gap = 18.1 mm)	35.9 keV
energy range (wiggler mode)	5.0 - 200.0 keV
energy range 1st harmonic (undulator mode)	0.2 - 4.0 keV
energy range 1st - 5th harmonics (undulator mode)	0.2 - 15.0 keV
on-axis peak brilliance at 2.7 keV (undulator mode)	$2.0 \times 10^{18}$ ph/sec/mrad <sup>2</sup> /mm <sup>2</sup> /0.1%bw
on-axis peak brilliance at 29.9 keV (wiggler mode at minimum gap = 18.1 mm)	$1.1 \times 10^{17}$ ph/sec/mrad <sup>2</sup> /mm <sup>2</sup> /0.1%bw
on-axis peak angular flux at 29.9 keV (wiggler mode at minimum gap = 18.1 mm)	$5.4 \times 10^{15}$ ph/sec/mrad <sup>2</sup> /0.1%bw
source size at critical energy $\Sigma_x$ $\Sigma_y$	359 $\mu$ m 21 $\mu$ m
source divergence at critical energy $\Sigma_{x'}$ (FWHM 1.2 mrad, non-Gaussian) $\Sigma_{y'}$	510 $\mu$ rad 47 $\mu$ rad

## Optics & Optical Performance

- slits
- power filters
- CARS-design Si(111) monochromator  
28.013 m from source (on orbit along  
0 mrad line)  
double bounce  
6.5–18.5 energy range w/ Si(111) crystal sets  
10<sup>-4</sup> energy resolution ( $\Delta E/E$ ) at 10 keV  
1.5 in. offset (fixed-exit, down bounce)  
water or liquid-nitrogen cooling

- CARS design cylindrical focusing mirror  
29.870 m from source (on orbit along  
0 mrad line)  
Si substrate, Rh coating  
water cooling at midplane  
4.1 mrad design angle  
horizontal focus: sagittal cylindrical figure  
vertical focus: bender

## Experiment Stations

- 14-ID-A  
white beam first optics enclosure  
10.5 m x 1.7 m x 2.8 m (L x W x H)
- 14-ID-B  
pink or monochromatic beam station  
4.5 m x 2.6 m x 2.8 m (L x W x H)  
MAD phasing  
microcrystallography  
virus/protein  
Laue  
time-resolved crystallography

## Detectors

- ADSC Q4, MAR345 and off-line image plate detectors
- NaI scintillation detectors
- Ge fluorescence detector

## Beamline Controls and Data Acquisition

- beamline controls: HP9000/778 HPUX 10.2 running EPICS via VME
- experiment: SGI 02 IRIX 6.3 running ADSC or Mar control software via in-house (kbscan>) software
- analysis: Fast SGI 02 IRIX 6.5 running DENZO, MOSFILM, DPS, etc.

## Beamline Support Equipment/Facilities

- cryo-coolers: Oxford CryoJet and CARS liquid-nitrogen/liquid-helium cooler
- collimators, filters, slits, beam-stop, CCD alignment camera
- beam position monitors (1  $\mu$ m resolution)
- beam flux monitors
- BL3 facility, sample-prep areas, cold room

## Beamline contacts:

Keith Moffat  
Wilfried Schildkamp  
Gary Navrotsky

tel 773-702-9950  
tel 630.252.0445  
tel 630.252.0444

moffat@cars.uchicago.edu  
schildkamp@cars.uchicago.edu  
navrotsky@cars.uchicago.edu

# Beamline 15-ID / ChemMatCARS-CAT

**Scientific focus:** Dynamic and structural condensed matter chemistry

**Scientific programs:** Time-dependent chemical crystallography scattering, anomalous scattering, micro crystallography; static and time-dependent surface scattering, dynamic protein diffraction, and small beam probing of complex structural polymers; interfacial and bulk studies using small- and wide-angle x-ray scattering; and coherent x-ray scattering of polymer/metal nanocomposites, nano-colloidal, and opaque materials.

## Insertion Device Source Characteristics (nominal)

source	Undulator A
period	3.30 cm
length	2.47 m
effective $K_{\max}$ (at minimum gap = 10.5 mm)	2.78
energy range 1st harmonic	2.9 - 13.0 keV
energy range 1st - 5th harmonics	2.9 - 45.0 keV
on-axis peak brilliance at 6.5 keV	$9.6 \times 10^{18}$ ph/sec/mrad <sup>2</sup> /mm <sup>2</sup> /0.1% bw
source size at 8.0 keV $\Sigma_x$ $\Sigma_y$	359 $\mu\text{m}$ 21 $\mu\text{m}$
source divergence at 8.0 keV $\Sigma_{x'}$ $\Sigma_{y'}$	24 $\mu\text{rad}$ 6.9 $\mu\text{rad}$

## Optics & Optical Performance

- Kohzu Seiki monochromator HLD-3  
3.1–22.5 keV energy range (for Si(111)  
and 25-mm offset)  
5°–40° angular range  
25–35 mm offset  
cryo-cooled Si or water-cooled diamond modes
- Oxford/SESIO vertically focusing mirror (A)  
water cooling  
white beam compatible
- Oxford/SESIO water-cooled 2nd steering mirror (B)
- modes of operation:  
1) monochromator w/ or w/out mirror(s)  
2) white beam mirror operation
- high-energy-resolution monochromator  
 $\Delta E/E \sim 10^{-5}$   
used with Kohzu monochromator
- in-station optics  
steering crystal or multilayer for surface science  
microfocusing optics

## Experiment Stations

- 15-ID-A  
white beam first optics enclosure
- 15-ID-B  
white beam station  
crystallography
- 15-ID-C  
pink beam station  
surface science
- 15-ID-D  
SAXS/WAXS  
10 m L x 5.8 m H

## Detectors

- Fuji image plates
- Bruker model 1500/1000 CCD
- Princeton Scientific model LN/CCD-1024SF CCD
- NaI
- ionization chambers
- avalanche photodiodes

## Beamline Controls and Data Acquisition

- EPICS and SPEC in addition to IDL, Windows NT and Sun Workstations running channel access with tools such as MEDM
- electronics VM- and NIM-based

## Beamline Support Equipment/Facilities

- Huber 6-circle (15-ID-B & -C)
- liquid surface spectrometer (15-ID-C)
- monochromatic beam chopper for time-resolved studies  
open time  $\sim 2.4 \mu\text{sec}$   
frequency  $\sim 1.3 \text{ kHz}$   
attenuation  $\sim 2 \times 10^{-7}$  at 33 keV
- Nd:YAG laser ( $\lambda=355 \text{ nm}$ , 400  $\mu\text{J}$  at 1000 Hz)


## Beamline contacts:

David Cookson  
Timothy Graber  
Yifei Jaski  
P. James Viccaro

tel 630.252.0468  
tel 630.252.0462  
tel 630.252.0461  
tel 630.252.0464

cookson@cars.uchicago.edu  
graber@cars.uchicago.edu  
ruan@cars.uchicago.edu  
viccaro@cars.uchicago.edu



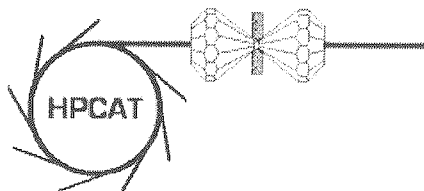

## High-Pressure CAT

---

### Overview

The High-Pressure Collaborative Access Team (HP-CAT) currently has three institutional members: the Carnegie Institution of Washington (managing member), Lawrence Livermore National Laboratory, and the University of Nevada, Las Vegas. These members pool resources and share expertise for the development, construction, and operation of sector 16. The purpose of HP-CAT is to provide its members with facilities for state-of-the-art educational and research activities and to assist the members and general users in operating these facilities; the goal is to advance high-pressure research with synchrotron radiation. The ultimate aim of HP-CAT is to set up a world-leading high-pressure research center at the APS and to make it accessible to the scientific community.

### Research Focus

Extreme pressure conditions radically alter the properties of materials. With the recent advances and developments in the technology of pressure vessels, which can now reach several hundreds of gigapascal (1 GPa = 10 kbar), the pressure variable promises to add a whole new dimension in fundamental physics, chemistry, earth and planetary sciences, and material science and technology. HP-CAT is dedicated to the scientific exploration of these fields with a common focus on the high-pressure dimension.

In ultrahigh pressure research, the power of an integrated approach is far greater than the sum of individual techniques. Comprehensive understanding of high-pressure phenomena relies on the combination of complementary measurements. HP-CAT integrates x-ray diffraction and spectroscopy with samples at high pressures and variable temperatures (from cryogenic conditions to thousands of degrees). These techniques include single-crystal structure refinement using energy-dispersive x-ray diffraction, high-resolution angle-dispersive polycrystalline x-ray diffraction, high-resolution x-ray emission spectroscopy, x-ray absorption spectroscopy, nuclear resonance forward scattering, inelastic x-ray scattering, and Compton scattering. The structural, vibrational, electronic, and magnetic properties of materials can thus be measured *in situ* for identifying novel phenomena and transitions at high pressures.

---

CAT contacts:	Ho-Kwang (David) Mao, CAT Dir.	tel 202.686.4454	mao@gl.ciw.edu
	Malcolm Nicol, HP-CAT Member	tel 702.895.1725	nicol@physics.unlv.edu
	Choong-Shik Yoo, HP-CAT Member	tel 925.422.5848	yoo1@llnl.gov
	Daniel Häusermann, Proj. Man.	tel 630.252.0427	daniel.hausermann@anl.gov
	Markus Schworer-Bohning, Beamline Scientist	tel 630.252.0411	schwore@aps.anl.gov

# Beamline 16-BM / HP-CAT

**Scientific focus:** Study of materials at high pressures in the fields of physics, chemistry, materials science, and planetary science

**Scientific programs:** Studies of materials under high pressure and variable temperature using diamond-anvil cells and other high-pressure devices by energy-dispersive diffraction from single crystals and powders and EXAFS

## Bending Magnet Source Characteristics (nominal)

source	APS bending magnet
critical energy	19.51 keV
energy range	1.0 - 100.0 keV
on-axis peak brilliance at 16.3 keV	$2.9 \times 10^{15}$ ph/sec/mrad <sup>2</sup> /mm <sup>2</sup> /0.1%bw
on-axis peak angular flux at 16.3 keV	$9.6 \times 10^{13}$ ph/sec/mrad <sup>2</sup> /0.1%bw
on-axis peak horizontal angular flux at 5.6 keV	$1.6 \times 10^{13}$ ph/sec/mradh/0.1%bw
source size at critical energy $\sum_x$ $\sum_y$	145 $\mu$ m 36 $\mu$ m
source divergence at critical energy $\sum_{x'}$ $\sum_{y'}$	6 mrad 47 $\mu$ rad

NOTE: The design parameters of beamline 16-BM are currently undergoing detailed reviews. Hence, the general characteristics given below will be refined and altered in the near future. Therefore, these parameters should be regarded as indicative only.

## Optics & Optical Performance

- water-cooled collimating mirror
- double-crystal monochromator  
water-cooled 1st crystal  
sagittally focusing 2nd crystal
- vertically focusing mirror
- small Kirkpatrick-Baez microfocusing mirrors

## Experiment Stations

- 16-BM-B  
white beam diffraction techniques using microfocusing
- 16-BM-C  
EXAFS experiments

## Detectors

- Ge and Si(Li) detectors
- other detectors to be defined

## Beamline Controls and Data Acquisition

- Windows NT-UNIX-Linux workstations running EPICS with VME
- SPEC, IDL, and other suitable control and data acquisition programs

## Beamline Support Equipment/Facilities

- custom-designed diffractometers and other sample handling stages
- diamond-anvil cells
- large-volume cell
- pressure measurement instrumentation
- cryostat
- optical spectrometers

**Beamline contacts:** Daniel Häusermann  
Markus Schwoerer-Bohning

tel 630.252.0427  
tel 630.252.0411

daniel.hausermann@anl.gov  
schwoere@aps.anl.gov



# Beamline 16-ID / HP-CAT

**Scientific focus:** Study of materials at high pressures in the fields of physics, chemistry, materials science, and planetary science

**Scientific programs:** Studies of materials under high pressure and variable temperature using diamond-anvil cells and other high-pressure devices by microfocus angle-dispersive diffraction from single crystals and powders, nuclear resonance forward scattering, inelastic scattering, Compton scattering, and emission and absorption spectroscopy

## Insertion Device Source Characteristics (nominal)

source	Undulator A
period	3.30 cm
length	2.47 m
effective $K_{\max}$ (at minimum gap = 10.5 mm)	2.78
energy range 1st harmonic	2.9 - 13.0 keV
energy range 1st - 5th harmonics	2.9 - 45.0 keV
on-axis peak brilliance at 6.5 keV	$9.6 \times 10^{18}$ ph/sec/mrad <sup>2</sup> /mm <sup>2</sup> /0.1% bw
source size at 8.0 keV $\sum_x$ $\sum_y$	359 $\mu\text{m}$ 21 $\mu\text{m}$
source divergence at 8.0 keV $\sum_{x'}$ $\sum_{y'}$	24 $\mu\text{rad}$ 6.9 $\mu\text{rad}$

NOTE: The design parameters of beamline 16-ID are currently undergoing detailed reviews. Hence, the general characteristics given below will be refined and altered in the near future. Therefore, these parameters should be regarded as indicative only.

## Optics & Optical Performance

- fixed-exit double-crystal monochromator  
diamond(111)  
water-cooled
- two-crystal branching monochromator  
diamond(111)  
water-cooled
- small Kirkpatrick-Baez microfocusing mirrors
- four-reflections / high-resolution monochromator
- vertically and horizontally focusing long mirrors

## Experiment Stations

- 16-ID-B  
monochromatic side-branch station for microbeam diffraction
- 16-ID-D  
white and monochromatic station for diffraction, spectroscopy, new techniques, and developments
- 16-ID-E  
white and monochromatic station for spectroscopy with high-resolution inelastic x-ray scattering spectrometer

## Detectors

- avalanche photodiode system for nuclear scattering
- Si PIN diode system for inelastic scattering
- area detector for angle-dispersive diffraction (to be defined)
- other detectors to be defined

## Beamline Controls and Data Acquisition

- Windows NT-UNIX-Linux workstations running EPICS with VME
- SPEC, IDL, and other suitable control and data acquisition programs

## Beamline Support Equipment/Facilities


- custom-designed diffractometers and other sample handling stages
- diamond-anvil cells
- large-volume cell
- pressure measurement instrumentation
- cryostat
- laser-heating system
- optical spectrometers

**Beamline contacts:** Daniel Häusermann  
Markus Schwoerer-Bohning

tel 630.252.0427  
tel 630.252.0411

daniel.hausermann@anl.gov  
schwoere@aps.anl.gov



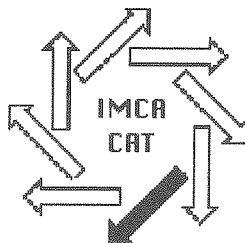

## Industrial Macromolecular Crystallography Association CAT

---

### Overview

The Industrial Macromolecular Crystallography Association (IMCA) is a consortium of crystallographic groups from 12 companies with major pharmaceutical research laboratories in the United States. IMCA-CAT comprises both IMCA and the Center for Synchrotron Radiation Research (CSRRI) at the Illinois Institute of Technology. The 12 corporate members of IMCA are Abbott Laboratories; Bristol-Myers Squibb Pharmaceutical Research Institute; Glaxo Wellcome Research Institute; Eli Lilly and Company; Merck and Company, Incorporated; Monsanto / Searle; Pfizer Pharmaceuticals; Pharmacia and Upjohn, Incorporated; The Procter and Gamble Company; Schering-Plough Research Institute; SmithKline Beecham Pharmaceuticals; and Three-Dimensional Pharmaceuticals. IMCA-CAT is unusual in that a large fraction of research conducted by CAT members is expected to be proprietary. Member scientists function as competitors in the discovery research that arises from their data-collection efforts; however, they collaborate in methods development and collectively fund the operation of the beamlines.

### Research Focus

Macromolecular crystallography is the primary research focus for IMCA, with the ultimate objective of facilitating the process of rational drug design. Traditional drug design methods involve testing thousands of compounds to determine biological activity. Rational drug design research attempts to elucidate biological processes in detail, which depends heavily on detailed knowledge of the atomic structure of the molecules involved. IMCA-CAT's insertion-device beamline is designed primarily for monochromatic and multiwavelength anomalous diffraction (MAD) experiments on small (<80 micrometer) protein crystals over an energy range of 6 to 20 keV. The bending-magnet line is designed mostly for monochromatic and MAD experiments on somewhat larger protein crystals over a similar energy range.

---

<b>CAT contacts:</b>	Andy Howard, Director	tel 630.252.0534	howard@iit.edu
	John Chrzas, Assoc. Director	tel 630.252.0522	john@metis.imca.aps.anl.gov
	Virginia Brown, CAT Coordinator	tel 630.252.0520	vbrown@metis.imca.aps.anl.gov

# Beamline 17-BM / IMCA-CAT

**Scientific focus:** Pharmaceutical macromolecular crystallography

**Scientific programs:** Structures of protein–ligand complexes, *de novo* protein structures, drug design, protein engineering, and crystallographic methods development

## Bending Magnet Source Characteristics (nominal)

source	APS bending magnet
critical energy	19.51 keV
energy range	1.0 - 100.0 keV
on-axis peak brilliance at 16.3 keV	$2.9 \times 10^{15}$ ph/sec/mrad <sup>2</sup> /mm <sup>2</sup> /0.1%bw
on-axis peak angular flux at 16.3 keV	$9.6 \times 10^{13}$ ph/sec/mrad <sup>2</sup> /0.1%bw
on-axis peak horizontal angular flux at 5.6 keV	$1.6 \times 10^8$ ph/sec/mradh/0.1%bw
source size at critical energy $\sum x$ $\sum y$	145 $\mu\text{m}$ 36 $\mu\text{m}$
source divergence at critical energy $\sum x'$ $\sum y'$	6 mrad 47 $\mu\text{rad}$

## Optics & Optical Performance

- Daresbury double-crystal constant off-set monochromator
  - 28 m from source
  - 6–20 keV energy range Si(111)
  - $10^{-4}$  energy resolution ( $\Delta E/E$ ) at 10 keV
  - 35 mm offset below orbital plane
  - water cooling
- sagittally bent 2nd monochromator crystal
- Daresbury vertically focusing mirror

## Experiment Stations

- 17-BM-A
  - white beam first optics enclosure
- 17-BM-B
  - white/monochromatic beam station
  - monochromatic macromolecular crystallography

## Detectors

- MAR345 image plate for crystallography
- fluorescence detector

## Beamline Controls and Data Acquisition

- controls: Sun and Linux systems running EPICS with VME “MX” software (locally developed), running on UNIX
- data acquisition: proprietary software from Mar

## Beamline Support Equipment/Facilities

- 4° chill room in wet lab
- user-accessible computers for data processing
- Oxford Instruments cryojet for sample cooling

**Beamline contacts:** Andrew Howard  
John Chrzas

tel 630.252.0534  
tel 630.252.0522

howard@iit.edu  
john@metis.imca.aps.anl.gov

# Beamline 17-ID / IMCA-CAT

**Scientific focus:** Pharmaceutical macromolecular crystallography

**Scientific programs:** Structures of protein- ligand complexes, *de novo* protein structures, drug design, protein engineering, and crystallographic methods development

## Insertion Device Source Characteristics (nominal)

source	Undulator A
period	3.30 cm
length	2.47 m
effective $K_{\max}$ (at minimum gap = 10.5 mm)	2.78
energy range 1st harmonic	2.9 - 13.0 keV
energy range 1st - 5th harmonics	2.9 - 45.0 keV
on-axis peak brilliance at 6.5 keV	$9.6 \times 10^{18}$ ph/sec/mrad <sup>2</sup> /mm <sup>2</sup> /0.1% bw
source size at 8.0 keV $\sum_x$ $\sum_y$	359 $\mu\text{m}$ 21 $\mu\text{m}$
source divergence at 8.0 keV $\sum_{x'}$ $\sum_{y'}$	24 $\mu\text{rad}$ 6.9 $\mu\text{rad}$

## Optics & Optical Performance

- Daresbury double-crystal constant off-set monochromator  
28 m from source  
6–20 keV energy range Si(111)  
 $10^{-4}$  energy resolution ( $\Delta E/E$ ) at 10 keV  
35 mm offset below orbital plane  
liquid-nitrogen cooling
- vertically focusing mirror under development

## Experiment Stations

- 17-ID-A  
white beam first optics enclosure
- 17-ID-B  
monochromatic beam station  
macromolecular crystallography  
3.5 m x 7 m

## Detectors

- Mar 165-mm single-element CCD
- Bicron fluorescence detector

## Beamline Controls and Data Acquisition

- controls: Sun and Linux systems running EPICS with VME “MX” software (locally developed), running on UNIX
- data acquisition: proprietary software from Mar

## Beamline Support Equipment/Facilities


- 4° chill room in wet lab
- Oxford cryo-stream system
- user-accessible computers for data processing

**Beamline contacts:** Andrew Howard  
John Chrzas

tel 630.252.0534  
tel 630.252.0522

howard@iit.edu  
john@metis.imca.aps.anl.gov



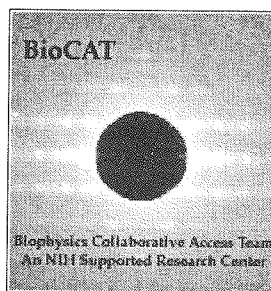

## Biophysics CAT

---

### Overview

The Biophysics CAT (Bio-CAT) is organized as a national research resource to be used primarily to study the structure of partially ordered biological molecules, complexes of biomolecules, and cellular structures under conditions similar to those present in living cells. The goal of this research is to determine the detailed mechanism of action of biological systems at the molecular level. The techniques used are x-ray fiber diffraction, x-ray scattering, x-ray absorption/emission spectroscopy, and advanced imaging techniques such as diffraction enhanced imaging. The research conducted at Bio-CAT's beamlines complements that done by other APS CATs (specifically, those that emphasize protein crystallography: BioCARS, IMCA-CAT, and SBC-CAT).

### Research Focus

Many biological macromolecules (e.g., enzymes) and complexes (e.g., viruses) can be crystallized and subjected to classic crystallographic analysis. However, in living cells, most biological systems are noncrystalline, and many biological structures cannot be studied in this manner because they cannot be crystallized, are transiently formed, or have structures or dynamical behaviors that change upon crystallization. For example, comparison of x-ray absorption spectra from crystalline and solution enzyme samples has shown that the processes of crystallization can change the microstructures sufficiently that crystallographic structural analysis alone would lead to incorrect interpretations of the enzymes' mechanisms of action.

All biological systems have some degree of spatial or dynamic order; this order can be probed by noncrystallographic diffraction and x-ray absorption. Of key interest to Bio-CAT scientists are the structures of small ordered domains such as single muscle fibers and connective tissues and studies of the time dependence of structural changes that occur as biomolecular complexes carry out their functions. Key research areas involve biological processes and systems of direct biomedical importance such as enzymes (particularly metalloproteins), DNA-binding proteins, proteins involved in gene expression, cell membranes, nerve cells, immune system components, the processes of cell transport, and cell motility.

---

<b>CAT contacts:</b>	Grant Bunker, Director	tel 312.567.3385	bunker@biocat1.iit.edu
	Thomas Irving, Deputy Dir.	tel 312.567.3489	irving@biocat1.iit.edu
	Clareen Krolik, CAT Coord.	tel 630.252.0549	krolik@bio.aps.anl.gov

# Beamline 18-ID / Bio-CAT

**Scientific focus:** Static and dynamic studies of partially ordered biological systems

**Scientific programs:** X-ray fiber diffraction, x-ray scattering, x-ray absorption/emission spectroscopy, and x-ray imaging and tomography

## Insertion Device Source Characteristics (nominal)

source	Undulator A
period	3.30 cm
length	2.47 m
effective $K_{\max}$ (at minimum gap = 10.5 mm)	2.78
energy range 1st harmonic	2.9 - 13.0 keV
energy range 1st - 5th harmonics	2.9 - 45.0 keV
on-axis peak brilliance at 6.5 keV	$9.6 \times 10^{18}$ ph/sec/mrad <sup>2</sup> /mm <sup>2</sup> /0.1% bw
source size at 8.0 keV $\sum_x$ $\sum_y$	359 $\mu\text{m}$ 21 $\mu\text{m}$
source divergence at 8.0 keV $\sum_{x'}$ $\sum_{y'}$	24 $\mu\text{rad}$ 6.9 $\mu\text{rad}$

## Optics & Optical Performance

- 3.5–35 keV energy range
- vertical beam size: 2 mm unfocused–  
0.02 mm focused
- horizontal beam size: 5 mm unfocused–  
0.01 mm focused
- Rosenbaum-Rock high-resolution monochromator #1  
6.5°–38° Bragg angle range  
1st crystal, Si(111), liquid-nitrogen cooling  
2nd crystal, Si(111), 10 mm wide (stabilized at 25°C, sagittally bent)  
7:1 demagnification  
35 mm beam offset (nominal)  
motorized tune, twist, & roll control  
(4 mm range, 50 nm resolution)

- Rosenbaum-Rock high-flux monochromator #2  
6.5°–38° Bragg angle range  
1st crystal, Si(400), liquid-nitrogen cooling  
2nd crystal, Si(400), 10 mm wide (stabilized at 25°C, sagittally bent)  
7:1 demagnification  
35 mm beam offset (nominal)  
motorized tune, twist, & roll control  
(4 mm range, 50 nm resolution)
- Rosenbaum-Rock vertical focusing mirror  
plane mirror substrate:  
Zerodur  
1020 mm x 95 mm x 38 mm  
2 Å rms roughness  
4  $\mu\text{rad}$  surface figure error, mounted  
Pt, none, and Pd coating stripes  
two motorized, encoded supports  
dynamic, independent bending mechanisms at both ends  
aberration correction via elliptical bending  
11:1 demagnification

## Experiment Station

- 18-ID-D  
monochromatic beam station  
12 m x 5 m x 3 m

## Detectors

- Fuji BAS 2500 image plate scanner
- 1024 x 1024 CCD detector with 60  $\mu\text{m}$  pixels
- fast-time-slicing plastic scintillator array
- Lytle fluorescence detector
- multilayer fluorescence analyzer

## Beamline Support Equipment/Facilities

- small-angle camera (up to 8 m camera length)
- Huber 4-circle (small), Huber 2-circle with quarter chi-segment
- Displex cryostat
- 3 ft x 5 ft optical table with five axes of motion

**Beamline contacts:** Grant Bunker  
Tom Irving  
Clareen Krolik

tel 312.567.3385  
tel 312.567.3489  
tel 630.252.0549

bunker@biocat1.iit.edu  
irving@biocat1.iit.edu  
krolik@bio.aps.anl.gov



## *Structural Biology Center*

### **Structural Biology Center CAT**

---

#### **Overview**

The Structural Biology Center CAT has constructed a national user facility at the APS for studies in macromolecular crystallography. Funded by the Department of Energy's Office of Biological and Environmental Research, the Structural Biology Center provides users with a full complement of instrumentation, ancillary facilities, software, and support staff to enhance the productivity of their research. The majority of beam time is allocated to Independent Investigators through a proposal process, with a small amount of time reserved for research and development activities by SBC-CAT staff members.

#### **Research Focus**

The principal focus of the SBC-CAT research program is macromolecular crystallography. Designed as a rapid-throughput facility, the Structural Biology Center provides researchers with the ability to use both standard crystallographic techniques and multiple energy anomalous dispersion (MAD) phasing methods.

One current major program focuses on the study of macromolecular assemblies, with an emphasis on molecular recognition and interaction between macromolecules. Specific projects include molecular chaperones (bacterial and archael chaperonins 60, human hsc70 chaperone), protein/nucleic acid complexes (trp repressor/operator systems), and multimeric enzymes.

---

<b>CAT contacts:</b>	Andrzej Joachimiak, Director WWW Site	tel 630.252.3926	andrzejj@anl.gov <a href="http://www.sbc.anl.gov">http://www.sbc.anl.gov</a>
----------------------	--	------------------	---

# Beamline 19-BM / SBC-CAT

**Scientific focus:** Macromolecular crystallography

**Scientific program:** Structural biology

## Bending Magnet Source Characteristics (nominal)

source	APS bending magnet
critical energy	19.51 keV
energy range	1.0 - 100.0 keV
on-axis peak brilliance at 16.3 keV	$2.9 \times 10^{15}$ ph/sec/mrad <sup>2</sup> /mm <sup>2</sup> /0.1%bw
on-axis peak angular flux at 16.3 keV	$9.6 \times 10^{13}$ ph/sec/mrad/0.1%bw
on-axis peak horizontal angular flux at 5.6 keV	$1.6 \times 10^{13}$ ph/sec/mradh/0.1%bw
source size at critical energy $\sum_x$ $\sum_y$	145 $\mu$ m 36 $\mu$ m
source divergence at critical energy $\sum_{x'}$ $\sum_{y'}$	6 mrad 47 $\mu$ rad

## Optics & Optical Performance

- 6–20 keV standard spectral range
- 0.11 mm hor. x 0.06 mm vert. FWHM focused spot size
- Rosenbaum-Rock double-crystal monochromator
  - water cooled
  - sagittally focusing 2nd crystal
  - 6:1 demagnification
  - 6°–40° Bragg angle range
  - 1st crystals: Si(111), Si(220), Si(331)
    - 25 mm wide ea., common cooler carrier
  - 2nd crystals: Si(111), Si(220), Si(331)
    - exchangeable, 25 mm x 125 mm
  - 35 mm beam offset (nominal)

- Rosenbaum-Rock vertical focusing mirror
  - 9:1 demagnification
  - Zerodur plane mirror substrate
  - 1020 mm x 100 mm x 38 mm
  - 2 Å rms roughness
  - 1  $\mu$ rad surface figure error
  - Pt, none, Pd coating stripes (35 mm wide ea.)
  - two motorized, encoded supports
  - dynamic, independent bending mechanism at both ends
  - aberration correction via elliptical bending

## Experiment Stations

- 19-BM-A
  - white beam first optics enclosure
- 19-BM-C
  - white beam optics enclosure
- 19-BM-D
  - monochromatic experiment station
  - kappa goniostat for macromolecular crystallography
  - guard slits
  - filter/shutter
  - detector support and positioner

## Detectors

- SBC1 3k x 3k CCD
  - built by ANL-ECT
  - 210 mm x 210 mm active area
  - 1.8 sec readout

## Beamline Controls and Data Acquisition

- Multiprocessor SG1 workstation, plus two UNIX work stations for data acquisition and data processing
- 3 HP workstations running EPICS, VME for beamline and detector control
- PMAC motor controller, software by ANL-ECT
- GUI for beamline control, data acquisition, and detector control by ANL-ECT

## Beamline Support Equipment/Facilities

- Rosenbaum-Rock miniaturized kappa goniostat
- high-magnification alignment cameras (two)
- Rosenbaum-Rock high-precision detector support and positioner
- liquid-nitrogen cryosystem sample cooler

**Beamline contacts:** Randy Alkire  
Stephan Ginell  
Andrzej Joachimiak

tel 630.252.3865      alkire@anl.gov  
tel 630.252.3972      ginell@anl.gov  
tel 630.252.3926      andrzejj@anl.gov

# Beamline 19-ID / SBC-CAT

**Scientific focus:** Macromolecular crystallography

**Scientific program:** Structural biology

## Insertion Device Source Characteristics (nominal)

source	Undulator A
period	3.30 cm
length	2.47 m
effective $K_{\max}$ (at minimum gap = 10.5 mm)	2.78
energy range 1st harmonic	2.9 - 13.0 keV
energy range 1st - 5th harmonics	2.9 - 45.0 keV
on-axis peak brilliance at 6.5 keV	$9.6 \times 10^{18}$ ph/sec/mrad <sup>2</sup> /mm <sup>2</sup> /0.1% bw
source size at 8.0 keV $\sum_x$ $\sum_y$	359 $\mu$ m 21 $\mu$ m
source divergence at 8.0 keV $\sum_{x'}$ $\sum_{y'}$	24 $\mu$ rad 6.9 $\mu$ rad

## Optics & Optical Performance

- 6–20 keV standard spectral range
- 0.08 mm hor. x 0.04 mm vert. FWHM focused beam size
- Rosenbaum-Rock monochromator #1
  - high-resolution double-crystal sagittal focusing
  - 1st crystal, Si(220), liquid-nitrogen cooled
  - 2nd crystal, Si(220), 10 mm wide (stabilized at 25°C; sagittally bent)
  - 7:1 demagnification
  - 6.5°–38° Bragg angle range
  - 35 mm beam offset (nominal)
  - motorized tune, twist, & roll control (4 mm range, 50 nm resolution)
- Rosenbaum-Rock monochromator #2
  - 1st crystal, Si(111), liquid-nitrogen cooled
  - 2nd crystal, Si(111), 10 mm wide, (stabilized at 25°C, sagittally bent)
  - 7:1 demagnification
  - 6.5°–38° Bragg angle range
  - 35 mm beam offset (nominal)
  - motorized tune, twist, & roll control (4 mm range, 50 nm resolution)

- Rosenbaum-Rock vertical focusing mirror
  - 11:1 demagnification
  - Zerodur plane mirror substrate
  - 1020 mm x 100 mm x 38 mm
  - 2 Å rms roughness
  - 1  $\mu$ rad surface figure error
  - Pt, none, and Pd coating stripes (35 mm wide ea.)
  - two motorized, encoded supports for angle and height adjustments
  - dynamic, independent bending mechanism at both ends
  - aberration correction via elliptical bending

## Experiment Stations

- 19-ID-A
  - white beam first optics enclosure
- 19-ID-C
  - white beam optics enclosure
- 19-ID-D
  - monochromatic experiment station
  - kappa goniostat for macromolecular crystallography
  - guard slits
  - filter/shutter
  - detector support and positioner

## Detectors

- SBC2 3k x 3k CCD
  - built by ANL-ECT
  - 210 mm x 210 mm active area
  - 1.8 sec readout

## Beamline Controls and Data Acquisition

- Multiprocessor SG1 workstation, plus two UNIX workstations for data acquisition
- Three HP workstations running EPICS, VME for beamline and detector control
- PMAC motor controller, software by ANL-ECT
- GUI for beamline control & data acquisition, detector control by ANL-ECT

## Beamline Support Equipment/Facilities

- Rosenbaum-Rock miniaturized kappa goniostat
- high-magnification alignment cameras (two)
- Rosenbaum-Rock high-precision detector support and positioner
- liquid-nitrogen cryosystem sample cooler

**Beamline contacts:** Randy Alkire  
Stephan Ginell  
Andrzej Joachimiak

tel 630.252.3865      alkire@anl.gov  
tel 630.252.3972      ginell@anl.gov  
tel 630.252.3926      andrzejj@anl.gov





## Pacific Northwest Consortium CAT

---

### Overview

The Pacific Northwestern Consortium CAT (PNC-CAT) includes a number of institutions in the Pacific Northwest regions of both the United States and Canada. Scientists from the lead institutions of the University of Washington, Pacific Northwest National Laboratories, and Simon Fraser University share several common research goals: specifically, a basic understanding of materials as related to their physical and chemical properties, the development of new strategies for environmental cleanup based on a fundamental understanding of the interaction of pollutants with the ambient environs, and structure-function properties of macromolecules.

### Research Focus

Several innovative instruments, including an analytical x-ray microscope with submicron resolution and capability to measure diffraction, x-ray absorption fine structure (XAFS), and diffraction anomalous fine structure (DAFS), will enable investigators to characterize heterogeneous materials in new ways. These measurements can be combined with two- and three-dimensional imaging using absorption and phase contrast. UHV and MBE facilities are available for *in situ* investigations using surface diffraction and XAFS. Capabilities for time-resolved measurements in the subnanosecond scale are planned. These techniques and instruments will be used to investigate materials science and environmental problems, on both a fundamental and applied level.

---

<b>CAT contacts:</b>	Edward Stern, Director	tel 206.543.2023	stern@phys.washington.edu
	Daryl Crozier, Associate Dir.	tel 604.291.4827	crozier@sfu.ca
	Steve Heald, Dir. of Construction	tel 630.252.9795	heald@aps.anl.gov

# Beamline 20-BM / PNC-CAT

**Scientific focus:** Materials science and environmental science

**Scientific programs:** Material science, environmental science, and surface science

## Bending Magnet Source Characteristics (nominal)

source	APS bending magnet
critical energy	19.51 keV
energy range	1.0 - 100.0 keV
on-axis peak brilliance at 16.3 keV	$2.9 \times 10^{15}$ ph/sec/mrad <sup>2</sup> /mm <sup>2</sup> /0.1%bw
on-axis peak angular flux at 16.3 keV	$9.6 \times 10^{13}$ ph/sec/mrad <sup>2</sup> /0.1%bw
on-axis peak horizontal angular flux at 5.6 keV	$1.6 \times 10^{13}$ ph/sec/mradh/0.1%bw
source size at critical energy $\sum_x$ $\sum_y$	145 $\mu$ m 36 $\mu$ m
source divergence at critical energy $\sum_{x'}$ $\sum_{y'}$	6 mrad 47 $\mu$ rad

## Optics & Optical Performance

- BESSRC-type monochromator
  - 29 m from source
  - 2–30 keV energy range
  - $1-2 \times 10^{-4}$  energy resolution ( $\Delta E/E$ )
  - 35 mm offset
  - water cooling
  - 2nd crystal : 4 mrad acceptance for sagittal focusing
- bent mirror (*planned*)
  - 26 m from source
  - 2–30 keV energy range
  - 2.7 mrad glancing angle

## Experiment Stations

- 20-BM-A  
white beam first optics enclosure
- 20-BM-B  
monochromatic general purpose station

## Detectors

- ionization chambers
- 13-element Ge
- CCDs (*planned*)

## Beamline Controls and Data Acquisition

- Sun, UNIX with EPICS/VME
- Windows NT with LabView
- SPEC (*planned*)

## Beamline Support Equipment/Facilities

- MBE/UHV surface chamber (*planned*)

**Beamline contacts:** Steve Heald  
Dale Brewé

tel 630.252.9795  
tel 630.252.0582

heald@aps.anl.gov  
brewé@pnc.aps.anl.gov

# Beamline 20-ID / PNC-CAT

**Scientific focus:** Materials science and environmental science

**Scientific programs:** Material science, environmental science, UHV/surface science-MBE growth, and x-ray microbeams

## Insertion Device Source Characteristics (nominal)

source	Undulator A
period	3.30 cm
length	2.47 m
effective $K_{\max}$ (at minimum gap = 10.5 mm)	2.78
energy range 1st harmonic	2.9 - 13.0 keV
energy range 1st - 5th harmonics	2.9 - 45.0 keV
on-axis peak brilliance at 6.5 keV	$9.6 \times 10^{18}$ ph/sec/mrad <sup>2</sup> /mm <sup>2</sup> /0.1% bw
source size at 8.0 keV $\sum_x$ $\sum_y$	359 $\mu$ m 21 $\mu$ m
source divergence at 8.0 keV $\sum_{x'}$ $\sum_{y'}$	24 $\mu$ rad 6.9 $\mu$ rad

## Optics & Optical Performance

- BESSRC high-heat-load monochromator
  - 31 m from source
  - 3–27 keV energy range, Si(111)
  - 6–50 keV energy range, Si(311)
  - $1-2 \times 10^{-4}$  energy resolution ( $\Delta E/E$ ), Si(111)
  - $0.5-1 \times 10^{-4}$  energy resolution ( $\Delta E/E$ ), Si(311)
  - 35 mm offset
  - liquid-nitrogen cooling
  - $\sim 10^{13}$  ph/sec flux at 10 keV, Si(111)
- toroidal focusing mirror
  - 34 m from source
  - 3–30 keV energy range
  - Al, Rh, Pt coatings (selectable)
  - 2.2–2.7 mrad glancing angle
  - $\sim 400 \times 200 \mu$ m (H,V) focal spot
- Kirkpatrick-Baez mirrors
  - 3–20 keV energy range
  - 1  $\mu$ m or 5  $\mu$ m focal spot (mirror dependent)
  - 4 mrad glancing angle
  - $\sim 10^{12}$  flux at 10 keV

- tapered capillaries
  - 3–30 keV energy range
  - $<0.5 \mu$ m focal spot
  - 150  $\mu$ m beam acceptance
  - $10^9-10^{10}$  ph/sec flux

## Experiment Stations

- 20-ID-A
  - white beam first optics enclosure
- 20-ID-B
  - monochromatic beam station
  - microbeams
  - XAFS
- 20-ID-C
  - monochromatic beam station
  - diffraction
  - surface science

## Detectors

- ionization chambers
- 13-element Ge detector
- NaI scintillation
- CCDs

## Beamline Controls and Data Acquisition

- Sun UNIX with EPICS/VME
- Windows NT with LabView
- SPEC

## Beamline Support Equipment/Facilities

- Huber surface diffractometer with UHV chamber (20-ID-C)
- MBE-capable surface XAFS, reflectivity and standing wave chamber (20-ID-B & C)
- diamond phase plate

**Beamline contacts:** Steve Heald  
Dale Brewé

tel 630.252.9795  
tel 630.252.0582

heald@aps.anl.gov  
brew@pnc.aps.anl.gov







## **Southeast Regional CAT**

---

### **Overview**

SER-CAT, which represents institutions from the states of Alabama, Florida, Georgia, Kentucky, Missouri, North Carolina, South Carolina, Tennessee, and Virginia, as well as the Intramural Program of the National Institutes of Health, plans to build and operate a sector at the APS. The SER-CAT mission is to provide researchers from its member institutions with timely access to tunable, high-intensity x-rays for x-ray diffraction structural studies of biomolecules. SER-CAT plans to provide its users with a complete research facility, including a full complement of instrumentation, software, and support staff for high-throughput data collection and processing. SER-CAT also plans to provide 25% of its beam time to outside users through its Independent Investigator program via a proposal process.

### **Research Focus**

The principal focus of the SER-CAT research program is macromolecular crystallography. The SER-CAT beamline designs are based on those of SBC-CAT; researchers will be able to use both traditional and multiple wavelength anomalous dispersion (MAD) techniques for structure determination and investigation.

The research interests of SER-CAT researchers are varied, ranging from DNA, protein and virus structures, to catalysts and basic materials science.

---

<b>CAT contacts:</b>	Bi-Cheng Wang, CAT Director	tel 706.542.1747	wang@BCL1.bmb.uga.edu
	Craig Smith, CAT Coordinator	tel 205.934.7233	smith@gold.cmc.uab.edu
	Gerold Rosenbaum, Project Dir.	tel 630.252.0643	rosenbaum@anl.gov
	John Unik, Project Administrator	tel 630.252.0640	junik@anl.gov

# Beamline 22-BM / SER-CAT

**Scientific focus:** Macromolecular crystallography

**Scientific program:** Structural biology

**Summary reflects projected performance parameters and planned equipment and hardware.**

## Bending Magnet Source Characteristics (nominal)

source	APS bending magnet
critical energy	19.51 keV
energy range	1.0 - 100.0 keV
on-axis peak brilliance at 16.3 keV	$2.9 \times 10^{15}$ ph/sec/mrad <sup>2</sup> /mm <sup>2</sup> /0.1%bw
on-axis peak angular flux at 16.3 keV	$9.6 \times 10^{13}$ ph/sec/mrad/0.1%bw
on-axis peak horizontal angular flux at 5.6 keV	$1.6 \times 10^{13}$ ph/sec/mradh/0.1%bw
source size at critical energy $\Sigma_x$ $\Sigma_y$	145 $\mu$ m 36 $\mu$ m
source divergence at critical energy $\Sigma_{x'}$ $\Sigma_{y'}$	6 mrad 47 $\mu$ rad

## Optics & Optical Performance

- 6–20 keV standard spectral range
- 0.11 mm hor. x 0.06 mm vert. FWHM focused spot size
- Rosenbaum-Rock double-crystal monochromator
  - water cooling
  - 1st crystals: Si(111), Si(220), Si(331)
    - 25 mm wide ea., common cooler carrier
  - 2nd crystals: Si(111), Si(220), Si(331)
    - sagital focusing
    - exchangeable
    - 125 mm L x 25 mm W
  - 6:1 demagnification
  - 6°–40° Bragg angle range
  - 35 mm beam offset (nominal)

- Rosenbaum-Rock vertical focusing mirror
  - 9:1 demagnification
  - ULE plane mirror substrate
  - 1020 mm x 100 mm x 38 mm
  - 2 Å rms roughness
  - 1  $\mu$ rad surface figure error
  - Pt, none, Pd coating stripes (35 mm wide ea.)
  - two motorized, encoded supports for angle and height adjustments
  - dynamic, independent bending mechanism at both ends
  - aberration correction via elliptical bending

## Experiment Stations

- 22-BM-A
  - white beam enclosure
- 22-BM-C
  - white beam optics enclosure
- 22-BM-D
  - monochromatic experiment station
  - kappa goniostat for macromolecular crystallography
  - adjustable collimator slits
  - filter/shutter
  - detector support and positioner

## Detectors

- CCD area detector

## Beamline Controls and Data Acquisition

- detector control and interface hardware and software to be defined
- beamline control and data acquisition software to be defined
- computers to be defined
- DC-servomotor; PMAC motor controller
- electrometer amplifiers; VME-based, computer-controlled V/F converter and scaler (ANL-ECT design)

## Beamline Support Equipment/Facilities

- Rosenbaum-Rock miniaturized kappa goniostat
- high-magnification alignment cameras (two)
- Rosenbaum-Rock high-precision detector support and positioner
- liquid-nitrogen cryosystem sample cooler

**Beamline contacts:** Gerd Rosenbaum  
John Gonczy

tel 630.252.0643  
tel 630.252.0642

rosenbaum@anl.gov  
gonczy@anl.gov

# Beamline 22-ID / SER-CAT

**Scientific focus:** Macromolecular crystallography

**Scientific program:** Structural biology

**Summary reflects projected performance parameters and planned equipment and hardware.**

## Insertion Device Source Characteristics (nominal)

source	Undulator A
period	3.30 cm
length	2.47 m
effective $K_{\max}$ (at minimum gap = 10.5 mm)	2.78
energy range 1st harmonic	2.9 - 13.0 keV
energy range 1st - 5th harmonics	2.9 - 45.0 keV
on-axis peak brilliance at 6.5 keV	$9.6 \times 10^{18}$ ph/sec/mrad <sup>2</sup> /mm <sup>2</sup> /0.1% bw
source size at 8.0 keV $\sum_x$ $\sum_y$	359 $\mu\text{m}$ 21 $\mu\text{m}$
source divergence at 8.0 keV $\sum_{x'}$ $\sum_{y'}$	24 $\mu\text{rad}$ 6.9 $\mu\text{rad}$

## Optics & Optical Performance

- 6–20 keV standard spectral range
- 0.08 mm hor. x 0.04 mm vert. FWHM focused beam size
- Rosenbaum-Rock monochromator
  - high-resolution double-crystal sagittal focusing
  - 1st crystal, Si(220), liquid-nitrogen cooling
  - 2nd crystal, Si(220), 10 mm wide (stabilized at 25°C, sagittally bent)
  - 7:1 demagnification
  - 6.5°–38° Bragg angle range
  - 35 mm beam offset (nominal)
  - motorized tune, twist, & roll control (4 mm range, 50 nm resolution)

- Rosenbaum-Rock vertical focusing mirror
  - 11:1 demagnification
  - ULE plane mirror substrate
  - 1020 mm x 100 mm x 38 mm
  - 2 Å rms roughness
  - 1  $\mu\text{rad}$  surface figure error
  - Pt, none, and Pd coating stripes (35 mm wide ea.)
  - two motorized, encoded supports
  - dynamic, independent bending mechanism at both ends
  - aberration correction via elliptical bending

## Experiment Stations

- 22-ID-A
  - white beam enclosure
- 22-ID-C
  - white beam optics enclosure
- 22-ID-D
  - monochromatic experiment station
  - kappa goniostat for macromolecular crystallography
  - adjustable collimator slits
  - filter/shutter
  - detector support and positioner

## Detectors

- CCD area detector

## Beamline Controls and Data Acquisition

- detector control and interface hardware and software to be defined
- beamline control and data acquisition software to be defined
- computers to be defined
- DC-servomotor; PMAC motor controller
- electrometer amplifiers; VME-based, computer-controlled V/F converter and scaler (ANL-ECT design)

## Beamline Support Equipment/Facilities

- Rosenbaum-Rock miniaturized kappa goniostat
- high-magnification alignment cameras (two)
- Rosenbaum-Rock high-precision detector support and positioner
- liquid-nitrogen cryosystem sample cooler

**Beamline contacts:** Gerd Rosenbaum  
John Goczy

tel 630.252.0643

tel 630.252.0642

rosenbaum@anl.gov

goczy@anl.gov



## SGX-CAT

### Structural GenomiX CAT

---

#### Overview

Structural GenomiX CAT (SGX-CAT) is designed for high-throughput crystallography in support of a proprietary research program in structural genomics. Funded entirely by Structural GenomiX (SGX), a San Diego-based biotechnology company, SGX-CAT is the first beamline facility at the Advanced Photon Source to be owned and operated by a single commercial entity. SGX is in the business of translating genome sequence information into protein structures for use in medicine, agriculture, and bioprocessing. The company will employ SGX-CAT to solve thousands of protein structures for use in drug and compound design by its customers. The SGX-CAT facility will include two undulator beamlines and a bending magnet beamline. The first undulator beamline is slated for completion in 2001.

#### Research Focus

The research focus of SGX-CAT is macromolecular crystallography. Structural GenomiX concentrates on proteins of interest to pharmaceutical, agricultural, and industrial enzyme customers, and it uses a structural genomics approach to achieve high throughput. Protein structures solved using SGX-CAT will be used to facilitate drug and compound design for Structural GenomiX customers. The SGX-CAT beamlines are designed to accommodate the high-throughput goals of Structural GenomiX. In particular, they will be optimized for rapid multiwavelength anomalous diffraction (MAD) experiments on small (~50 micron) selenomethionyl protein crystals.



## COM-CAT

### Commercial CAT

---

#### Overview

The Commercial Beamline Collaborative Access Team (COM-CAT) will be specializing in the application of synchrotron radiation to industrial problems. COM-CAT plans to operate a general-purpose beamline as an analytical service laboratory, providing both materials analysis and data interpretation for problems that require synchrotron radiation techniques for their solution. Funds for beamline construction were provided by the State of Illinois. Beamline operations will ultimately be funded on a fee-for-service basis.

#### Research Focus

Analyses based on synchrotron radiation are generally applicable to problems in materials science. Industrial companies have used synchrotron facilities to examine a wide variety of materials, including polymers, catalysts, candidate drugs, and biomolecules. COM-CAT expects to provide a convenient means to use synchrotron analysis capabilities to industrial organizations who either do not have current access to a beamline or require rapid turnaround of samples.

COM-CAT will, in its initial stages, construct and operate a single undulator beamline. To provide as broad a range of analytical capabilities as possible, the experimental facilities will consist of three permanent stations devoted to spectroscopy, single-crystal diffraction, small-angle scattering, and x-ray fluorescence. Zone plates will be used when microfocusing is required. To maximize use of the available beam, each table will have a permanent basic configuration. To use a downstream station, the upstream tables will be lowered below the flight path of the x-ray beam.

COM-CAT is also involved in outreach to the general industrial community. It has organized a workshop on industrial applications of synchrotron radiation and will provide on-site introductions to synchrotron radiation for potential customers.

---

<b>CAT contacts:</b>	Kevin D'Amico, Director	tel 630.252.3959	kdamico@aps.anl.gov
	Stephen Wasserman, Deputy Dir.	tel 630.252.3527	srw@anl.gov
	Cindy Doran, Secretary	tel 630.252.5564	doran@aps.anl.gov

# Beamline 32-ID / COM-CAT

**Scientific focus:** Materials characterization

**Scientific programs:** Analytical materials characterization and capability for commercial customers

## Insertion Device Source Characteristics (nominal)

source	Undulator A
period	3.30 cm
length	2.47 m
effective $K_{\max}$ (at minimum gap = 10.5 mm)	2.78
energy range 1st harmonic	2.9 - 13.0 keV
energy range 1st - 5th harmonics	2.9 - 45.0 keV
on-axis peak brilliance at 6.5 keV	$9.6 \times 10^{18}$ ph/sec/mrad <sup>2</sup> /mm <sup>2</sup> /0.1% bw
source size at 8.0 keV $\sum_x$ $\sum_y$	359 $\mu\text{m}$ 21 $\mu\text{m}$
source divergence at 8.0 keV $\sum_{x'}$ $\sum_{y'}$	24 $\mu\text{rad}$ 6.9 $\mu\text{rad}$

## Optics & Optical Performance

- L5-90 standard component primary aperture
- Kohzu HLD-3 monochromator
- horizontally focusing mirror  
Pd, Pt coating stripes

## Experiment Stations

- 32-ID-A  
white beam first optics enclosure
- 32-ID-B  
white beam station  
monochromatic experiments  
spectroscopy  
scattering  
crystallography

## Detectors

- Lytle detector
- multi-element energy dispersive detector
- ionization chambers
- scintillation counters
- CCD area detector

## Beamline Controls and Data Acquisition

- Linux PC running SPEC

## Beamline Support Equipment/Facilities

- 6-circle diffractometer; spectroscopy instruments
- crystallography station

**Beamline contacts:** Kevin D'Amico  
Stephen Wasserman

tel 630.252.3959  
tel 630.252.3527

kdamico@aps.anl.gov  
srw@anl.gov





## A University-National Laboratory-Industry CAT

---

### Overview

The UNI-CAT is a collaboration between scientists from the University of Illinois, the Oak Ridge National Laboratory, the National Institute of Standards and Technology, and UOP Research and Development. The UNI-CAT mission is to instrument and operate x-ray research facilities at the APS that provide advanced x-ray techniques to a diverse scientific community. This goal is achieved through the development of two sectors at the APS that are designed for cutting-edge x-ray studies in materials, physics, chemistry, biology, and geology. This is a multipurpose scattering facility capable of high-resolution scattering with excellent energy resolution and beam-focusing optics. Special capabilities in the UNI-CAT sectors also include a dedicated surface/interface diffraction station equipped with a molecular-beam epitaxy and chemical-vapor deposition facility, an ultra-small-angle x-ray scattering apparatus, an x-ray microscope, a topography station, a dedicated insertion-device beamline for microfocus and coherent x-ray diffraction, instrumentation for time-resolved x-ray scattering experiments, and reactive environment x-ray absorption spectroscopy.

### Research Focus

The primary research areas of UNI-CAT members encompass materials sciences and condensed matter physics, and include techniques such as structural crystallography, diffuse x-ray scattering, magnetic x-ray scattering, ultra-small-angle x-ray scattering, x-ray microscopy, millivolt resolution spectroscopy, surface and interface scattering, absorption spectroscopy, x-ray topography, microbeam techniques, coherent x-ray diffraction, and time-resolved techniques. These research tools permit UNI-CAT scientists and collaborators to explore fundamental structure/property relationships in bulk solids, at surfaces, and at internal interface boundaries.

---

<b>CAT contacts:</b>	Haydn Chen, Director	tel 217.244.4666	unicat@uimrl7.mrl.uiuc.edu
	Paul Zschack, Associate Dir.	tel 630.252.0860	zschack@anl.gov

# Beamline 33-BM / UNI-CAT

**Scientific focus:** Materials science and condensed matter physics

**Scientific programs:** Materials Science, ceramic science, phase transitions, surface science, thin-film structure and growth, and materials physics

## Bending Magnet Source Characteristics (nominal)

source	APS bending magnet
critical energy	19.51 keV
energy range	1.0 - 100.0 keV
on-axis peak brilliance at 16.3 keV	$2.9 \times 10^{15}$ ph/sec/mrad <sup>2</sup> /mm <sup>2</sup> /0.1%bw
on-axis peak angular flux at 16.3 keV	$9.6 \times 10^{13}$ ph/sec/mrad <sup>2</sup> /0.1%bw
on-axis peak horizontal angular flux at 5.6 keV	$1.6 \times 10^{13}$ ph/sec/mradh/0.1%bw
source size at critical energy $\sum_x$ $\sum_y$	145 $\mu$ m 36 $\mu$ m
source divergence at critical energy $\sum_{x'}$ $\sum_{y'}$	6 mrad 47 $\mu$ rad

## Optics & Optical Performance

- 4–40 keV energy range
- collimating mirror
- double-crystal monochromator
- focusing mirror

## Experiment Stations

- 33-BM-B  
EXAFS  
topography station
- 33-BM-C  
general purpose scattering station

## Detectors

- NaI scintillation counters
- gas-filled proportional counters (Xe, Ar)
- ionization chambers

## Beamline Controls and Data Acquisition

- Sun UNIX running EPICS with VME, SPEC

## Beamline Support Equipment/Facilities

- topography station (33-BM-B)
- EXAFS table (33-BM-B)
- Huber 4-circle diffractometer (33-BM-C)

**Beamline contact:** Paul Zschack

tel 630.252.0860

zschack@anl.gov

# Beamline 33-ID / UNI-CAT

**Scientific focus:** Material Science and condensed matter physics

**Scientific programs:** Material science, ceramic science, phase transitions, surface science, thin-films structure and growth, and materials physics

## Insertion Device Source Characteristics (nominal)

source	Undulator A
period	3.30 cm
length	2.47 m
effective $K_{\max}$ (at minimum gap = 10.5 mm)	2.78
energy range 1st harmonic	2.9 - 13.0 keV
energy range 1st - 5th harmonics	2.9 - 45.0 keV
on-axis peak brilliance at 6.5 keV	$9.6 \times 10^{18}$ ph/sec/mrad <sup>2</sup> /mm <sup>2</sup> /0.1% bw
source size at 8.0 keV $\sum_x$ $\sum_y$	359 $\mu\text{m}$ 21 $\mu\text{m}$
source divergence at 8.0 keV $\sum_{x'}$ $\sum_{y'}$	24 $\mu\text{rad}$ 6.9 $\mu\text{rad}$

## Optics & Optical Performance

- focused beam size 280  $\mu\text{m}$  hor. x 100  $\mu\text{m}$  vert.
- PSL double-crystal monochromator
  - 45 m from source
  - 4.0–40.0 keV energy range Si(111)
  - $10^{-4}$  energy resolution ( $\Delta E/E$ )
  - 20–35 mm beam offset
  - liquid-nitrogen cooling
  - variable sagittal focus
- vertical focusing mirror
  - 2.0–4.0 mrad angle of incidence
  - Pt, Si, Rh coating stripes
  - variable vertical focus
- harmonic rejection mirror
  - 2.0–4.0 mrad angle of incidence
  - Pt, Si, Rh coating stripes

## Experiment Stations

- 33-ID-D
  - general purpose scattering
  - USAXS station
- 33-ID-E
  - surface and interface scattering station

## Detectors

- NaI scintillation counters
- ionization chambers
- gas-filled proportional counters (Xe, Ar)

## Beamline Controls and Data Acquisition

- Sun UNIX running EPICS with VME, SPEC

## Beamline Support Equipment/Facilities

- Huber 4-circle diffractometer (33-ID-D)
- Newport 6-circle goniometer w/ kappa geometry (33-ID-D)
- USAXS (33-ID-D)
- UHV surface/interface scattering chamber (33-ID-E)

**Beamline contact:** Paul Zschack

tel 630.252.0860

zschack@anl.gov

# Beamline 34-ID / UNI-CAT

**Scientific focus:** Microfocussing and coherent x-ray scattering

**Scientific programs:** Coherent x-ray diffraction, microbeam scattering, and x-ray microprobe

## Insertion Device Source Characteristics (nominal)

source	Undulator A
period	3.30 cm
length	2.47 m
effective $K_{\max}$ (at minimum gap = 10.5 mm)	2.78
energy range 1st harmonic	2.9 - 13.0 keV
energy range 1st - 5th harmonics	2.9 - 45.0 keV
on-axis peak brilliance at 6.5 keV	$9.6 \times 10^{18}$ ph/sec/mrad <sup>2</sup> /mm <sup>2</sup> /0.1% bw
source size at 8.0 keV $\sum_x$ $\sum_y$	359 $\mu\text{m}$ 21 $\mu\text{m}$
source divergence at 8.0 keV $\sum_{x'}$ $\sum_{y'}$	24 $\mu\text{rad}$ 6.9 $\mu\text{rad}$

## Optics & Optical Performance

- beamline under design
- beam-splitting mirrors
- monochromators

## Experiment Stations

- 34-ID-C  
coherent beam diffraction station
- 34-ID-E  
microbeam scattering station

## Detectors

- various CCD cameras
- NaI scintillation counters

## Beamline Controls and Data Acquisition

- To be determined

## Beamline Support Equipment/Facilities

- Kappa diffractometer with surface science capability (34-ID-C)
- microbeam scattering table (34-ID-E)

**Beamline contacts:** Paul Zschack  
Ian Robinson  
Gene Ice

tel 630.252.0860  
tel 217.244.2949  
tel 423.574.2744

zschack@anl.gov  
robinson@uimrl7.mrl.uiuc.edu  
icege@ornl.gov

# Directory

	BUILDING	SECTOR(S)	PHONE	FAX	E-MAIL
APS User Office	401	n/a	630.252.9090	630.252.9250	apsuser@aps.anl.gov
BESSRC-CAT	433	11, 12	630.252.1553	630.252.0365	cowan@anl.gov
Bio-CAT	435	18	630.252.0549	630.252.0545	krollk@bio.aps.anl.gov
CARS-CAT	434	13, 14, 15	630.252.0450	630.252.0443	franklin@cars.uchicago.edu
CMC-CAT	433	9	630.252.0320	630.252.0339	gog@anl.gov
COM-CAT	438	32	630.252.0840	630.252.0855	doran@aps.anl.gov
DND-CAT	432	5	630.252.0222	630.252.0226	jhopf@nwu.edu
HP-CAT	434	16 (under devel.)	630.252.0427	630.252.0443	daniel.hausemann@anl.gov
IMCA-CAT	435	17	630.252.0520	630.252.0521	vbrown@metis.imca.aps.anl.gov
IMMW-CAT	432	8	630.252.0281	630.252.0282	larry@lindy.mit.edu
MHATT-CAT	432	7	630.252.0209	630.252.0279	williams@mail.mhattach.howard.edu
MR-CAT	433	10	219.631.7219	219.631.5952	bunker@nd.edu
$\mu$ -CAT	432	6	630.252.0247	630.252.0240	drobinsn@iastate.edu
PNC-CAT	435	20	630.252.9795	630.252.0580	heald@aps.anl.gov
SBC-CAT	435	19	630.252.0560	630.252.0564	mfcner@anl.gov
SER-CAT	363A	22 (under devel.)	630.252.0640	630.252.0641	junik@anl.gov
SGX-CAT	438	31 (under devel.)	630.252.3959	630.252.0855	kdamico@aps.anl.gov
SRI-CAT	431	1, 2, 3, 4	630.252.0160	630.252.0161	ljs@aps.anl.gov
UNI-CAT	438	33, 34	630.252.0860	630.252.0862	zschack@anl.gov

To call an insertion-device beamline, dial 630.252.18XX

To call a bending-magnet beamline, dial 630.252-17XX

(where XX is the sector number)